

Продукция



Мочеточниковые стенты



Нефростомические наборы



Эпцистостомические наборы



Проводники



Мочеточниковые катетеры





Катетеры Фолея



Детская урология


Россия:

Санкт-Петербург
Новорощинская ул. д.4
Офис 714-1
office@inmed.spb.ru

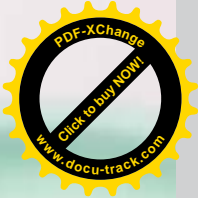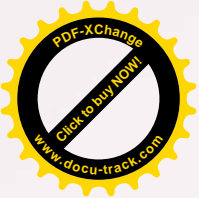

URETERAL STENTS



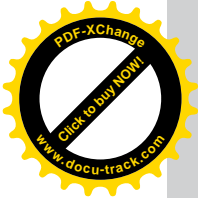

URETERAL STENTS

Variations overview

YELLOW-STAR

- Polyurethane (PUR)
- Phosphorylcholine (PC) coating to minimise the risk of infection and encrustation
- Very suitable for long term indwelling time



WHITE-STAR

- Polyurethane (PUR)
- Standard stent – available in many different versions



GREEN-STAR

- Polyurethane (PUR)
- Special soft material for high patient comfort



URETERAL STENTS

Introduction and general information

Introduction

On the following pages, we would like to provide you with a comprehensive portfolio of our ureteral stents. Regardless of the indication, you are sure to find the appropriate product here.

Material

Polyurethane (PUR)

All UROTECH ureteral stents are produced from aliphatic polyurethane (PUR) that has been specifically developed for medical applications. Unlike other materials, PUR shows minimal adhesion to blood cells and has a very low-grade inflammatory effect on tissue due to its physical and chemical properties. PUR becomes soft in a moist environment at body temperature, significantly limiting the risks of

traumatic injury to surrounding tissue.

Experimental tests and practical experience have shown that PUR is one of the least thrombogenic and most tissue-compatible materials for implantable medical products.

Indwelling time

All UROTECH ureteral stents are approved for long-term use and can remain inside a patient for several months. The additional layer of phosphorylcholine (PC) on the ureteral stents significantly reduces the risk of encrustation and infection.

However, due to the heterogeneity of indications and individual differences between patients, regular follow-up examinations should be performed by the treating physician.

Length measurements

In accordance with international norms (ASTM F-1828), the length of our ureteral stents is always indicated without including the J-tip – this equals the shaft length (Fig. 1).


Variations


In addition to our large catalog of standard ureteral stents, you can order many different custom variations based on your clinical needs.

E.g. stents:

- ◆ with attached suture
- ◆ without drainage openings in shaft

Fig. 1


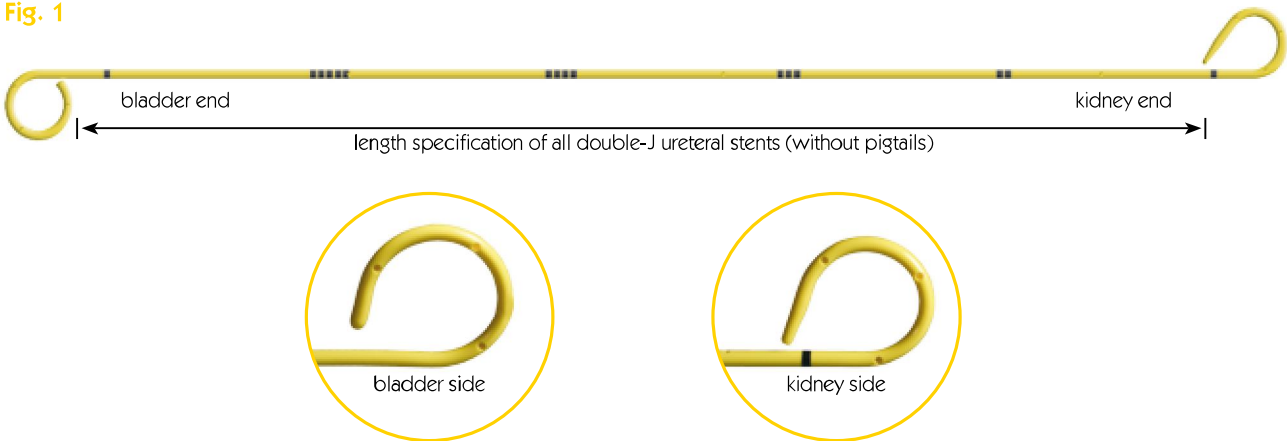


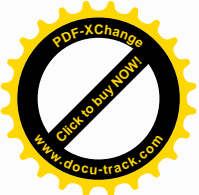

URETERAL STENTS

YELLOW-STAR




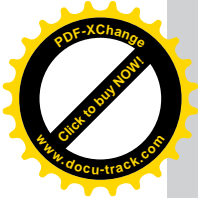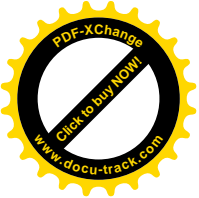

YELLOW-STAR Double-J Ureteral Stents-Set

PC-coated

SET COMPONENTS

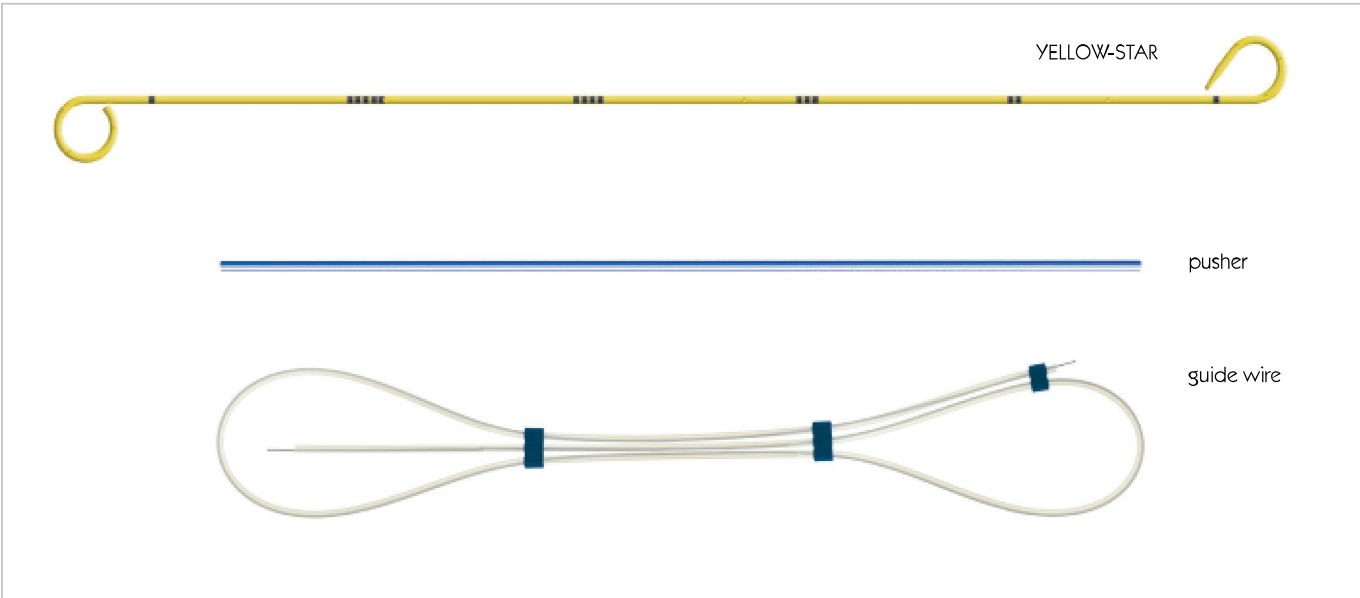

Ureteral Stent

- ◆ Polyurethane (PUR)
- ◆ PC coating to minimise risk of infection and encrustation, resulting in extended indwelling time
- ◆ Drainage openings in pigtails and shaft
- ◆ Open or closed tip
- ◆ Markings for exact positioning
- ◆ Radiopaque

Pusher

- ◆ Length: 40 cm



Guide Wire

- ◆ PTFE-coated

Clamp

- ◆ Only in sets with closed tip



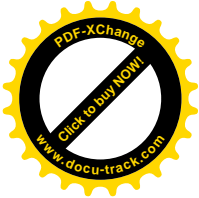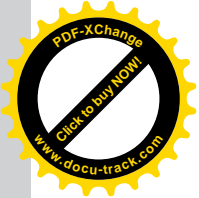



URETERAL STENTS

YELLOW-STAR Double-J Ureteral Stents-Set

PC-coated

Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24	EP-414724	EP-364724	125	0.028"	straight	standard
4.7	26	EP-414726	EP-364726	125	0.028"	straight	standard
06	24	EP-410624	EP-360624	125	0.035"	straight	standard
06	26	EP-410626	EP-360626	125	0.035"	straight	standard
06	28	EP-410628	EP-360628	125	0.035"	straight	standard
06	30	EP-410630	EP-360630	125	0.035"	straight	standard
07	24	EP-410724	EP-360724	125	0.035"	straight	standard
07	26	EP-410726	EP-360726	125	0.035"	straight	standard
07	28	EP-410728	EP-360728	125	0.035"	straight	standard
07	30	EP-410730	EP-360730	125	0.035"	straight	standard
08	24	EP-410824	EP-360824	125	0.035"	straight	standard
08	26	EP-410826	EP-360826	125	0.035"	straight	standard
08	28	EP-410828	EP-360828	125	0.035"	straight	standard
08	30	EP-410830	EP-360830	125	0.035"	straight	standard


Closed tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24	EP-404724	EP-394724	125	0.028"	straight	standard
4.7	26	EP-404726	EP-394726	125	0.028"	straight	standard
06	24	EP-400624	EP-390624	125	0.035"	straight	standard
06	26	EP-400626	EP-390626	125	0.035"	straight	standard
06	28	EP-400628	EP-390628	125	0.035"	straight	standard
06	30	EP-400630	EP-390630	125	0.035"	straight	standard
07	24	EP-400724	EP-390724	125	0.035"	straight	standard
07	26	EP-400726	EP-390726	125	0.035"	straight	standard
07	28	EP-400728	EP-390728	125	0.035"	straight	standard
07	30	EP-400730	EP-390730	125	0.035"	straight	standard
08	24	EP-400824	EP-390824	125	0.035"	straight	standard
08	26	EP-400826	EP-390826	125	0.035"	straight	standard
08	28	EP-400828	EP-390828	125	0.035"	straight	standard
08	30	EP-400830	EP-390830	125	0.035"	straight	standard

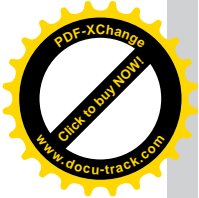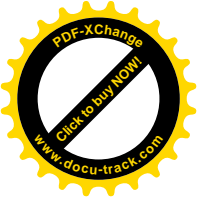

YELLOW-STAR Double-J Ureteral Stent-Set

PC-coated, flexible length

SET COMPONENTS

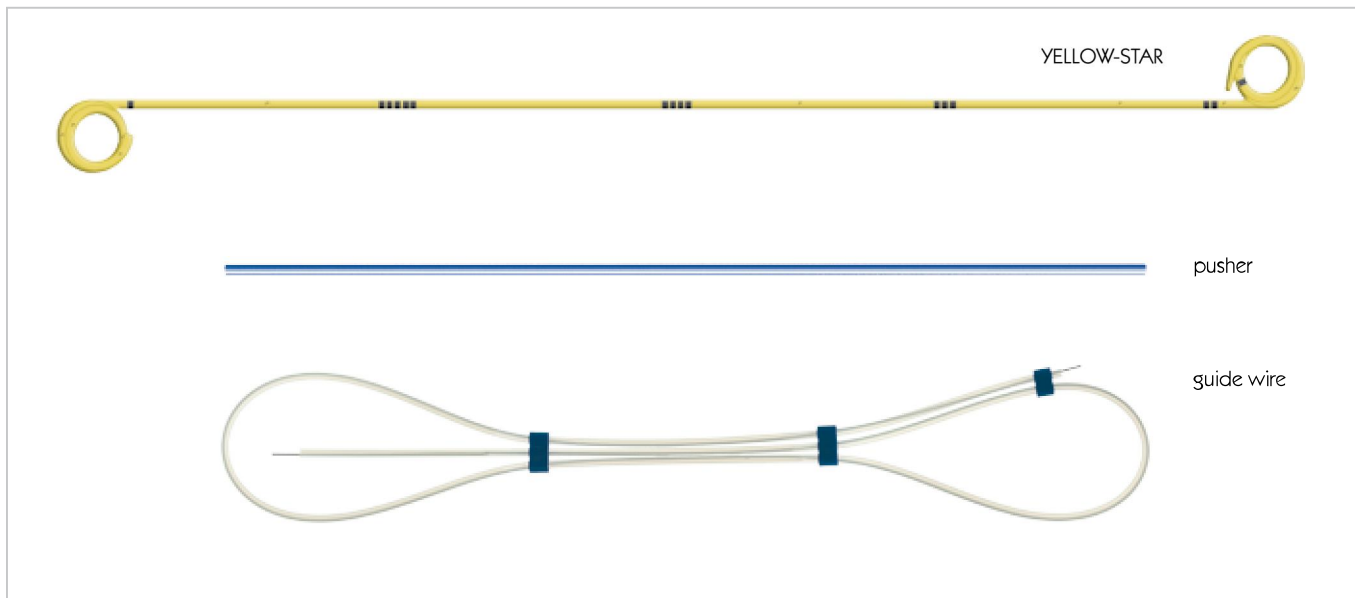


Ureteral Stent

- ♦ Polyurethane (PUR)
- ♦ PC coating to minimise risk of infection and encrustation, resulting in extended indwelling time
- ♦ Multiple coils on both ends for flexible shaft length from 22 to 30 cm
- ♦ Less risk of dislocation
- ♦ Drainage openings in pigtails and shaft
- ♦ Open or closed tip
- ♦ Markings for exact positioning
- ♦ Radiopaque



bladder side

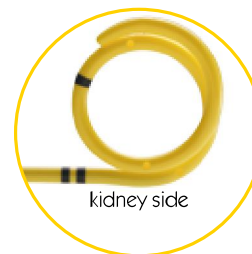

kidney side

Pusher

- ♦ Length: 40 cm


Guide Wire

- ♦ PTFE-coated

Clamp

- ♦ Only in sets with closed tip



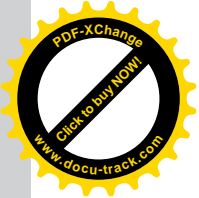


URETERAL STENTS

YELLOW-STAR Double-J Ureteral Stent-Set PC-coated, flexible length

Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	22-30	EP-414722-30-BS	EP-364722-30-BS	125	0.028"	straight	standard
06	22-30	EP-410622-30-BS	EP-360622-30-BS	125	0.035"	straight	standard
07	22-30	EP-410722-30-BS	EP-360722-30-BS	125	0.035"	straight	standard
08	22-30	EP-410822-30-BS	EP-360822-30-BS	125	0.035"	straight	standard


Closed tip							
Size (FR/CH)	Length (cm)		Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	22-30		EP-394722-30-BS	125	0.028"	straight	standard
06	22-30		EP-390622-30-BS	125	0.035"	straight	standard
07	22-30		EP-390722-30-BS	125	0.035"	straight	standard
08	22-30		EP-390822-30-BS	125	0.035"	straight	standard

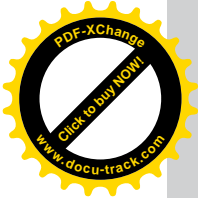

YELLOW-STAR Double-J Ureteral Stents-Set

PC-coated, steerable

SET COMPONENTS

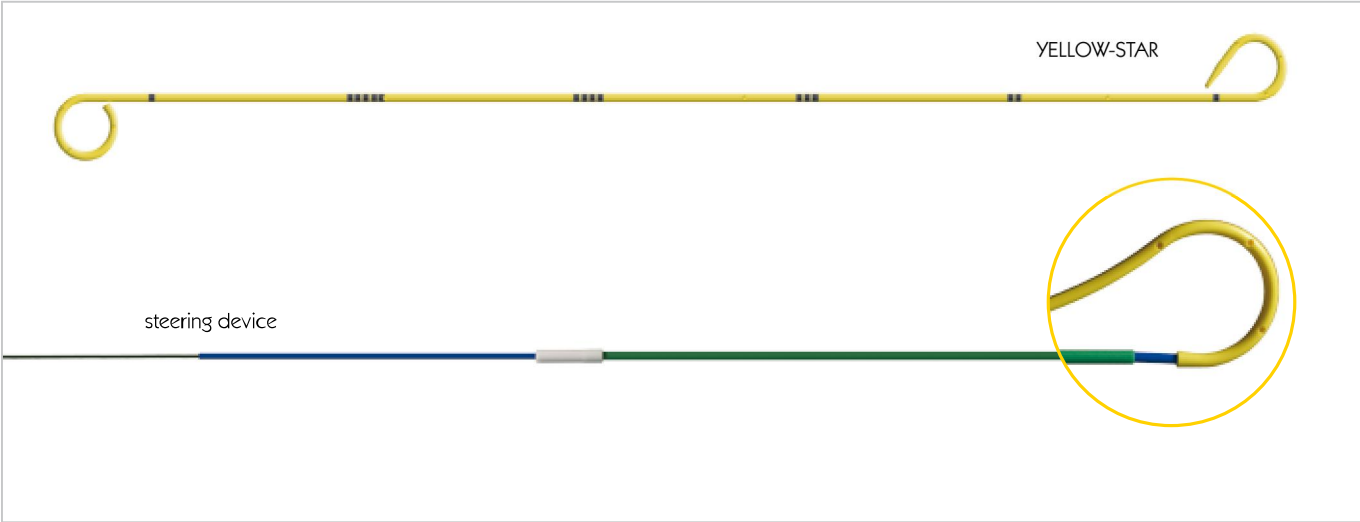

Ureteral Stent

- ◆ Polyurethane (PUR)
- ◆ PC coating to minimise risk of infection and encrustation, resulting in extended indwelling time
- ◆ Drainage openings in pigtails and shaft
- ◆ Open or closed tip
- ◆ Markings for exact positioning
- ◆ Radiopaque

Steering Device

- ◆ Pusher
- ◆ Disconnecter
- ◆ Length: 53 cm



Guide Wire

- ◆ PTFE-coated

Clamp

- ◆ Only in sets with closed tip



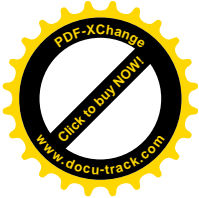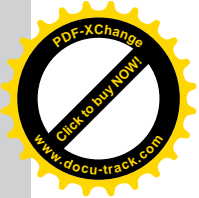


YELLOW-STAR Double-J Ureteral Stents-Set

PC-coated, steerable

Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24	EP-474724	EP-484724	145	0.028"	straight	standard
4.7	26	EP-474726	EP-484726	145	0.028"	straight	standard
06	24	EP-470624	EP-480624	150	0.035"	straight	standard / movable core
06	26	EP-470626	EP-480626	150	0.035"	straight	standard / movable core
06	28	EP-470628	EP-480628	150	0.035"	straight	standard / movable core
06	30	EP-470630	EP-480630	150	0.035"	straight	standard / movable core
07	24	EP-470724	EP-480724	150	0.035"	straight	standard / movable core
07	26	EP-470726	EP-480726	150	0.035"	straight	standard / movable core
07	28	EP-470728	EP-480728	150	0.035"	straight	standard / movable core
07	30	EP-470730	EP-480730	150	0.035"	straight	standard / movable core
08	24	EP-470824	EP-480824	150	0.035"	straight	standard / movable core
08	26	EP-470826	EP-480826	150	0.035"	straight	standard / movable core
08	28	EP-470828	EP-480828	150	0.035"	straight	standard / movable core
08	30	EP-470830	EP-480830	150	0.035"	straight	standard / movable core


Closed tip							
Size (FR/CH)	Length (cm)		Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24		EP-494724	125	0.028"	straight	standard
4.7	26		EP-494726	125	0.028"	straight	standard
06	24		EP-490624	125	0.035"	J-tip	standard / movable core
06	26		EP-490626	125	0.035"	J-tip	standard / movable core
06	28		EP-490628	125	0.035"	J-tip	standard / movable core
06	30		EP-490630	125	0.035"	J-tip	standard / movable core
07	24		EP-490724	125	0.035"	J-tip	standard / movable core
07	26		EP-490726	125	0.035"	J-tip	standard / movable core
07	28		EP-490728	125	0.035"	J-tip	standard / movable core
07	30		EP-490730	125	0.035"	J-tip	standard / movable core
08	24		EP-490824	125	0.035"	J-tip	standard / movable core
08	26		EP-490826	125	0.035"	J-tip	standard / movable core
08	28		EP-490828	125	0.035"	J-tip	standard / movable core
08	30		EP-490830	125	0.035"	J-tip	standard / movable core

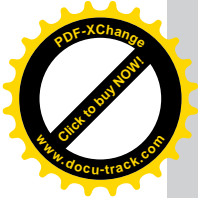

YELLOW-STAR Double-J Ureteral Stent-Set

PC-coated, steerable with coupling

SET COMPONENTS

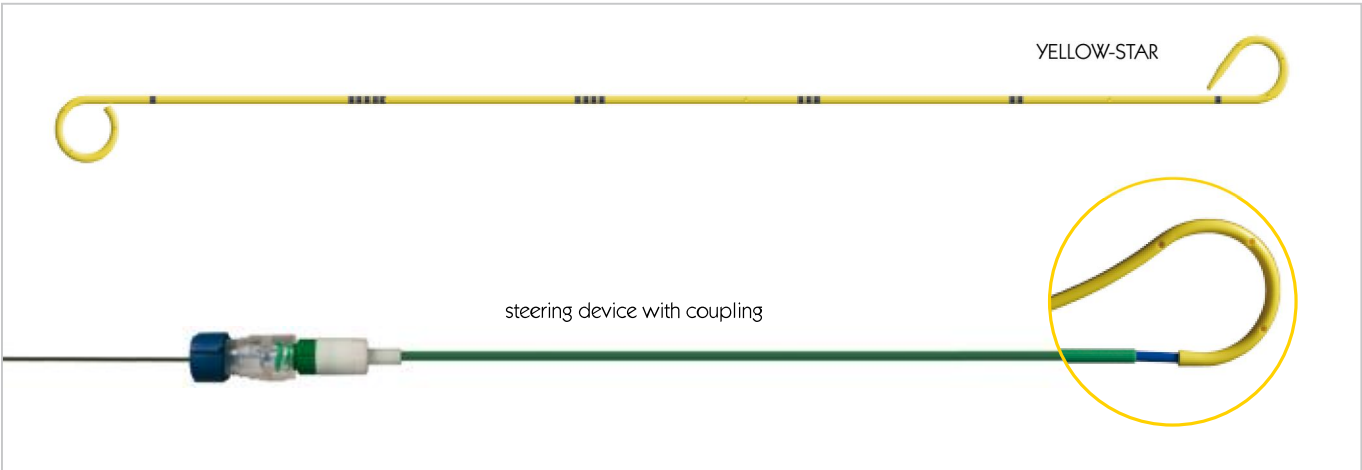

Ureteral Stent



- ♦ Polyurethane (PUR)
- ♦ PC coating to minimise risk of infection and encrustation, resulting in extended indwelling time
- ♦ Drainage openings in pigtails and shaft
- ♦ Open or closed tip
- ♦ Markings for exact positioning
- ♦ Radiopaque

Steering Device with Coupling

- ♦ Pusher
- ♦ Disconnecter
- ♦ Coupling to fasten guide wire
- ♦ Length: 43 cm

Guide Wire

- ♦ PTFE-coated

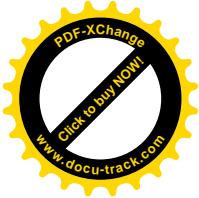
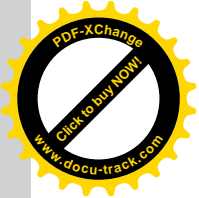


YELLOW-STAR Double-J Ureteral Stent-Set

PC-coated, steerable with coupling


Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24	EP-574724	EP-584724	145	0.028"	straight	standard
4.7	26	EP-574726	EP-584726	145	0.028"	straight	standard
06	24	EP-570624	EP-580624	150	0.035"	straight	standard / movable core
06	26	EP-570626	EP-580626	150	0.035"	straight	standard / movable core
06	28	EP-570628	EP-580628	150	0.035"	straight	standard / movable core
06	30	EP-570630	EP-580630	150	0.035"	straight	standard / movable core
07	24	EP-570724	EP-580724	150	0.035"	straight	standard / movable core
07	26	EP-570726	EP-580726	150	0.035"	straight	standard / movable core
07	28	EP-570728	EP-580728	150	0.035"	straight	standard / movable core
07	30	EP-570730	EP-580730	150	0.035"	straight	standard / movable core
08	24	EP-570824	EP-580824	150	0.035"	straight	standard / movable core
08	26	EP-570826	EP-580826	150	0.035"	straight	standard / movable core
08	28	EP-570828	EP-580828	150	0.035"	straight	standard / movable core
08	30	EP-570830	EP-580830	150	0.035"	straight	standard / movable core



Closed tip							
Size (FR/CH)	Length (cm)		Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24		EP-594724	125	0.028"	straight	standard
4.7	26		EP-594726	125	0.028"	straight	standard
06	24		EP-590624	125	0.035"	J-tip	standard / movable core
06	26		EP-590626	125	0.035"	J-tip	standard / movable core
06	28		EP-590628	125	0.035"	J-tip	standard / movable core
06	30		EP-590630	125	0.035"	J-tip	standard / movable core
07	24		EP-590724	125	0.035"	J-tip	standard / movable core
07	26		EP-590726	125	0.035"	J-tip	standard / movable core
07	28		EP-590728	125	0.035"	J-tip	standard / movable core
07	30		EP-590730	125	0.035"	J-tip	standard / movable core
08	24		EP-590824	125	0.035"	J-tip	standard / movable core
08	26		EP-590826	125	0.035"	J-tip	standard / movable core
08	28		EP-590828	125	0.035"	J-tip	standard / movable core
08	30		EP-590830	125	0.035"	J-tip	standard / movable core

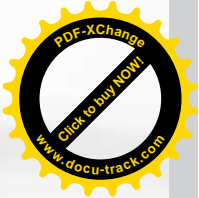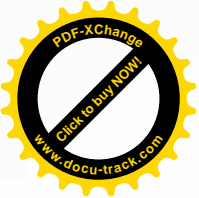

URETERAL STENTS


WHITE-STAR

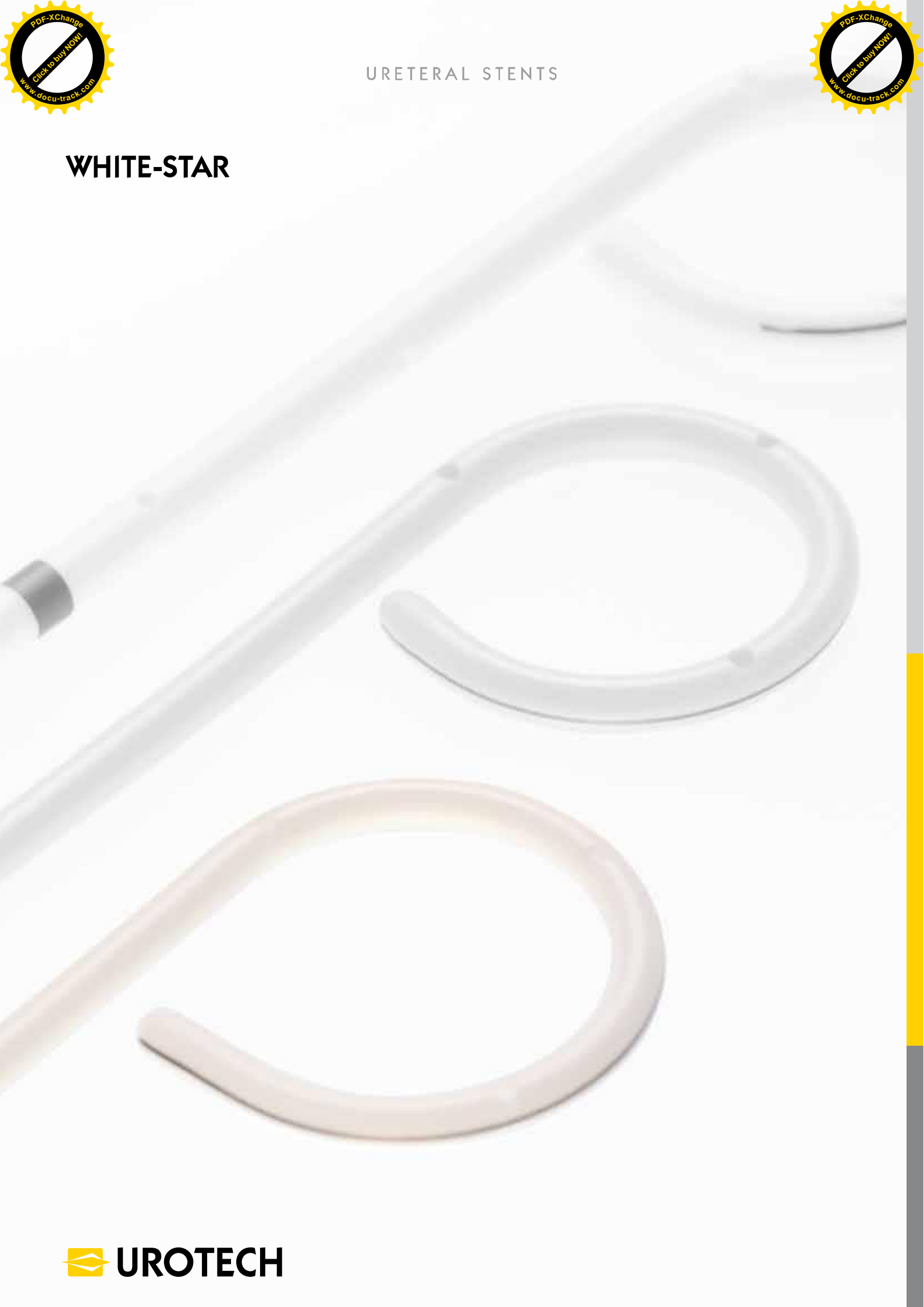

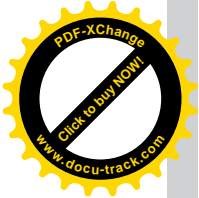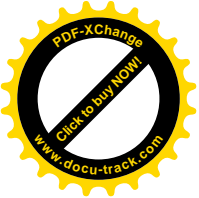

WHITE-STAR Double-J Ureteral Stents-Set

SET COMPONENTS

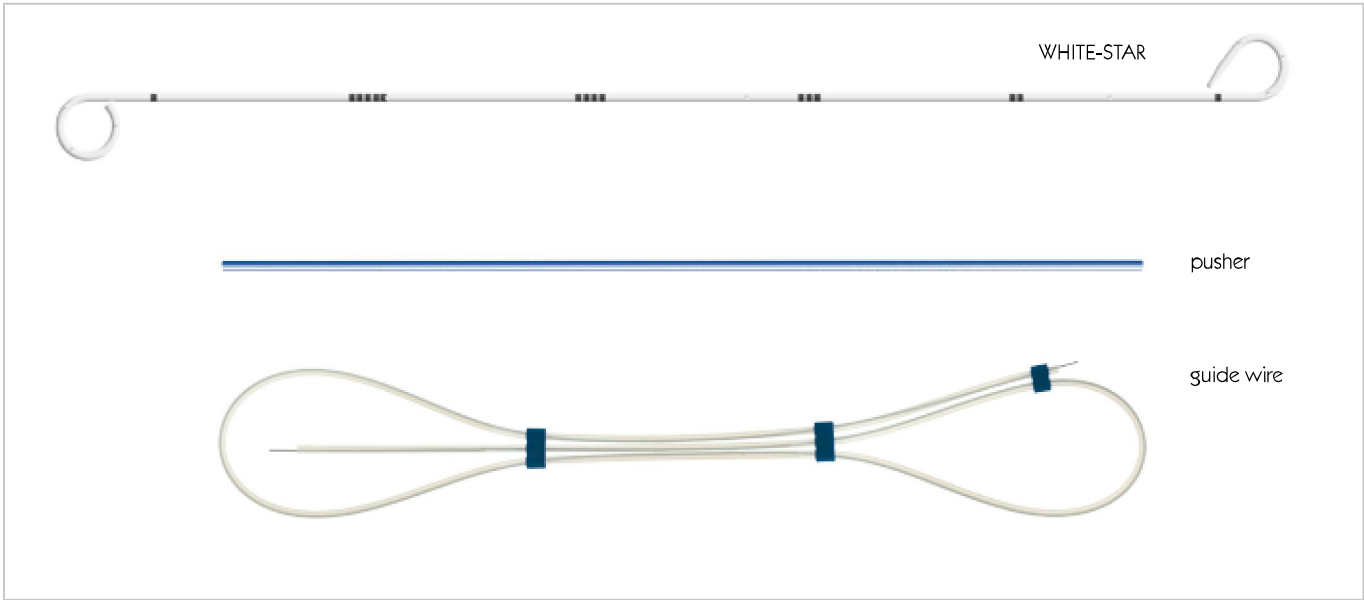

Ureteral Stent

- ◆ Polyurethane (PUR)
- ◆ Drainage openings in pigtails and shaft
- ◆ Open or closed tip
- ◆ Markings for exact positioning
- ◆ Radiopaque

Pusher

- ◆ Length: 40 cm



Guide Wire

- ◆ PTFE-coated

Clamp

- ◆ 2 clamps in sets with closed tip



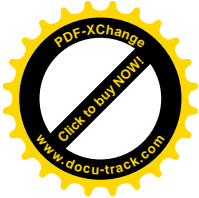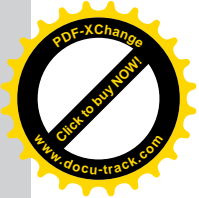


URETERAL STENTS

WHITE-STAR Double-J Ureteral Stents-Set

Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24	EN-410524	EN-360524	125	0.028"	straight	standard
4.7	26	EN-410526	EN-360526	125	0.028"	straight	standard
06	24	EN-410624	EN-360624	125	0.035"	straight	standard
06	26	EN-410626	EN-360626	125	0.035"	straight	standard
06	28	EN-410628	EN-360628	125	0.035"	straight	standard
06	30	EN-410630	EN-360630	125	0.035"	straight	standard
07	24	EN-410724	EN-360724	125	0.035"	straight	standard
07	26	EN-410726	EN-360726	125	0.035"	straight	standard
07	28	EN-410728	EN-360728	125	0.035"	straight	standard
07	30	EN-410730	EN-360730	125	0.035"	straight	standard
08	24	EN-410824	EN-360824	125	0.035"	straight	standard
08	26	EN-410826	EN-360826	125	0.035"	straight	standard
08	28	EN-410828	EN-360828	125	0.035"	straight	standard
08	30	EN-410830	EN-360830	125	0.035"	straight	standard


Closed tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24	EN-400524	EN-300524	125	0.028"	straight	standard
4.7	26	EN-400526	EN-300526	125	0.028"	straight	standard
06	24	EN-400624	EN-300624	125	0.035"	straight	standard
06	26	EN-400626	EN-300626	125	0.035"	straight	standard
06	28	EN-400628	EN-300628	125	0.035"	straight	standard
06	30	EN-400630	EN-300630	125	0.035"	straight	standard
07	24	EN-400724	EN-300724	125	0.035"	straight	standard
07	26	EN-400726	EN-300726	125	0.035"	straight	standard
07	28	EN-400728	EN-300728	125	0.035"	straight	standard
07	30	EN-400730	EN-300730	125	0.035"	straight	standard
08	24	EN-400824	EN-300824	125	0.035"	straight	standard
08	26	EN-400826	EN-300826	125	0.035"	straight	standard
08	28	EN-400828	EN-300828	125	0.035"	straight	standard
08	30	EN-400830	EN-300830	125	0.035"	straight	standard

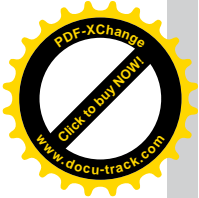

WHITE-STAR Double-J Ureteral Stents-Set

Flexible length

SET COMPONENTS

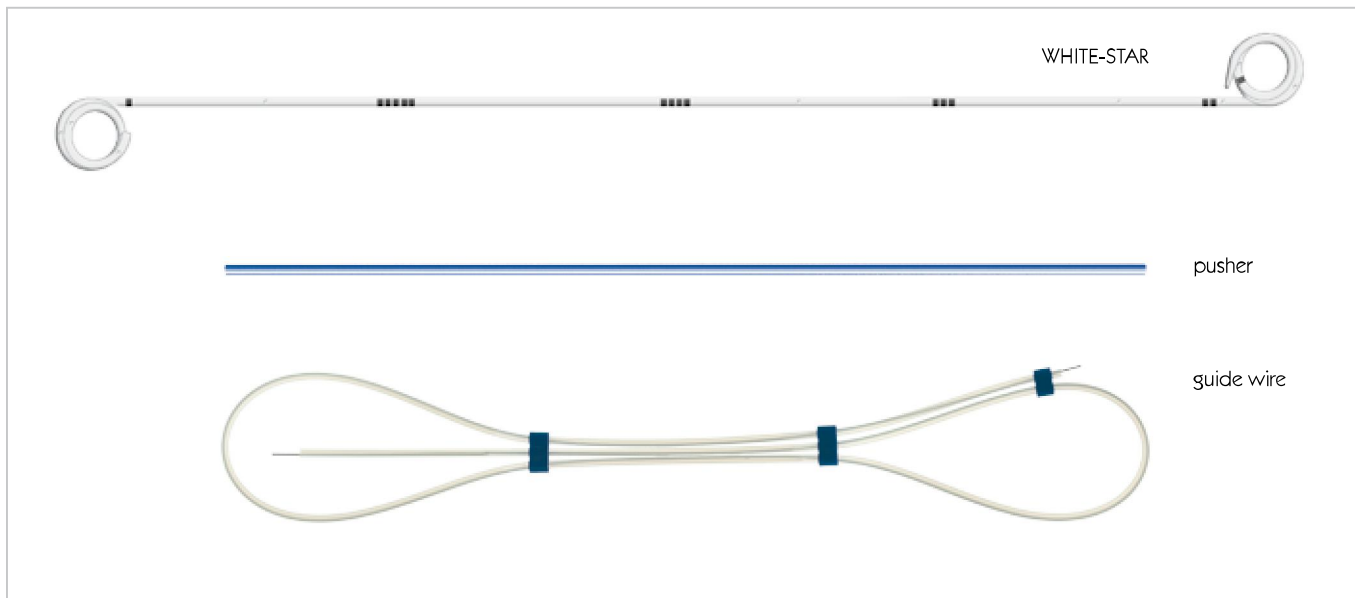

Ureteral Stent

- ♦ Polyurethane (PUR)
- ♦ Multiple coils on both ends for flexible shaft length from 22 to 30 cm
- ♦ Less risk of dislocation
- ♦ Drainage openings in pigtails and shaft
- ♦ Open or closed tip
- ♦ Markings for exact positioning
- ♦ Radiopaque



Pusher

- ♦ Length: 40 cm


Guide Wire

- ♦ PTFE-coated

Clamp

- ♦ 2 clamps in sets with closed tip



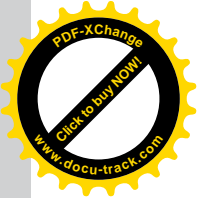


URETERAL STENTS

WHITE-STAR Double-J Ureteral Stents-Set Flexible length


Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	22-30	EN-410622-30-BS	EN-360622-30-BS	125	0.035"	straight	standard
07	22-30	EN-410722-30-BS	EN-360722-30-BS	125	0.035"	straight	standard
08	22-30	EN-410822-30-BS	EN-360822-30-BS	125	0.035"	straight	standard

Closed tip							
Size (FR/CH)	Length (cm)		Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	22-30		EN-390622-30-BS	125	0.035"	straight	standard
07	22-30		EN-390722-30-BS	125	0.035"	straight	standard
08	22-30		EN-390822-30-BS	125	0.035"	straight	standard

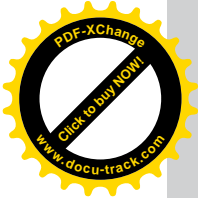

WHITE-STAR Double-J Ureteral Stents-Set Steerable

SET COMPONENTS

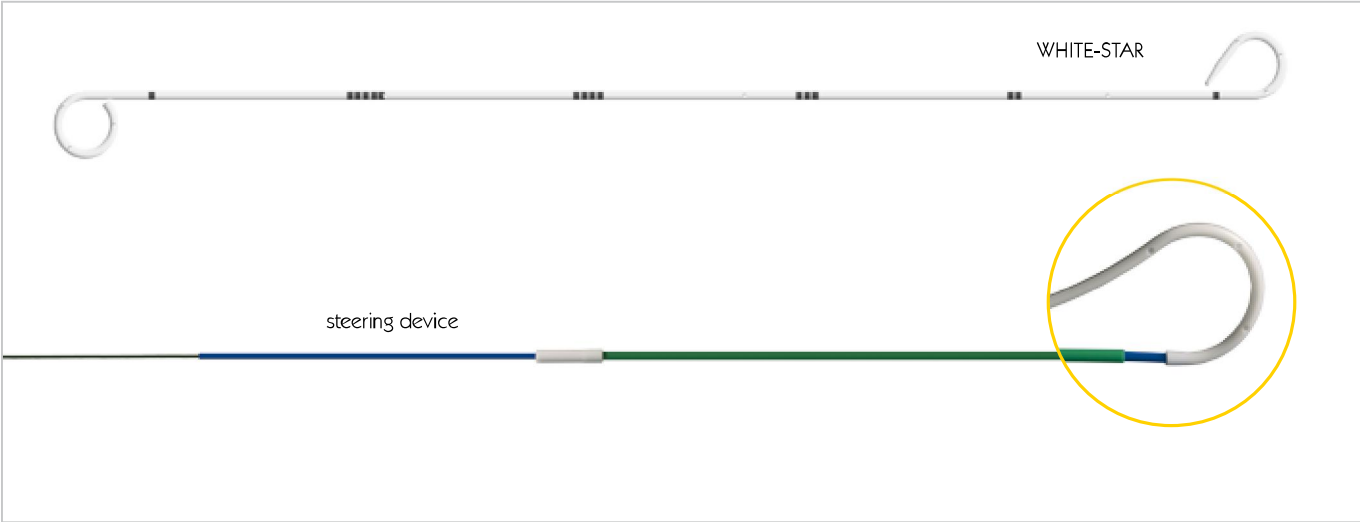

Ureteral Stent

- ♦ Polyurethane (PUR)
- ♦ Drainage openings in pigtails and shaft
- ♦ Open or closed tip
- ♦ Markings for exact positioning
- ♦ Radiopaque

Steering Device

- ♦ Pusher
- ♦ Disconnecter
- ♦ Length: 53 cm



Guide Wire

- ♦ PTFE-coated

Clamp

- ♦ 2 clamps in sets with closed tip



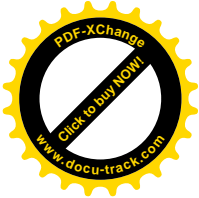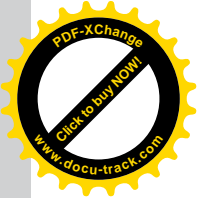



URETERAL STENTS

WHITE-STAR Double-J Ureteral Stents-Set

Steerable

Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24	ES-360524	ES-370524	145	0.028"	straight	standard
4.7	26	ES-360526	ES-370526	145	0.028"	straight	standard
06	24	ES-360624	ES-370624	150	0.035"	straight	standard / movable core
06	26	ES-360626	ES-370626	150	0.035"	straight	standard / movable core
06	28	ES-360628	ES-370628	150	0.035"	straight	standard / movable core
06	30	ES-360630	ES-370630	150	0.035"	straight	standard / movable core
07	24	ES-360724	ES-370724	150	0.035"	straight	standard / movable core
07	26	ES-360726	ES-370726	150	0.035"	straight	standard / movable core
07	28	ES-360728	ES-370728	150	0.035"	straight	standard / movable core
07	30	ES-360730	ES-370730	150	0.035"	straight	standard / movable core
08	24	ES-360824	ES-370824	150	0.035"	straight	standard / movable core
08	26	ES-360826	ES-370826	150	0.035"	straight	standard / movable core
08	28	ES-360828	ES-370828	150	0.035"	straight	standard / movable core
08	30	ES-360830	ES-370830	150	0.035"	straight	standard / movable core


Closed tip							
Size (FR/CH)	Length (cm)		Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24		ES-380524	125	0.028"	straight	standard
4.7	26		ES-380526	125	0.028"	straight	standard
06	24		ES-380624	125	0.035"	J-tip	standard / movable core
06	26		ES-380626	125	0.035"	J-tip	standard / movable core
06	28		ES-380628	125	0.035"	J-tip	standard / movable core
06	30		ES-380630	125	0.035"	J-tip	standard / movable core
07	24		ES-380724	125	0.035"	J-tip	standard / movable core
07	26		ES-380726	125	0.035"	J-tip	standard / movable core
07	28		ES-380728	125	0.035"	J-tip	standard / movable core
07	30		ES-380730	125	0.035"	J-tip	standard / movable core
08	24		ES-380824	125	0.035"	J-tip	standard / movable core
08	26		ES-380826	125	0.035"	J-tip	standard / movable core
08	28		ES-380828	125	0.035"	J-tip	standard / movable core
08	30		ES-380830	125	0.035"	J-tip	standard / movable core

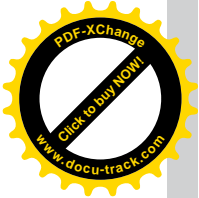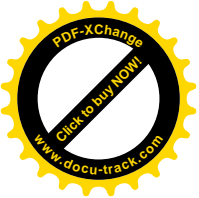

WHITE-STAR Double-J Ureteral Stent-Set

Steerable with coupling

SET COMPONENTS

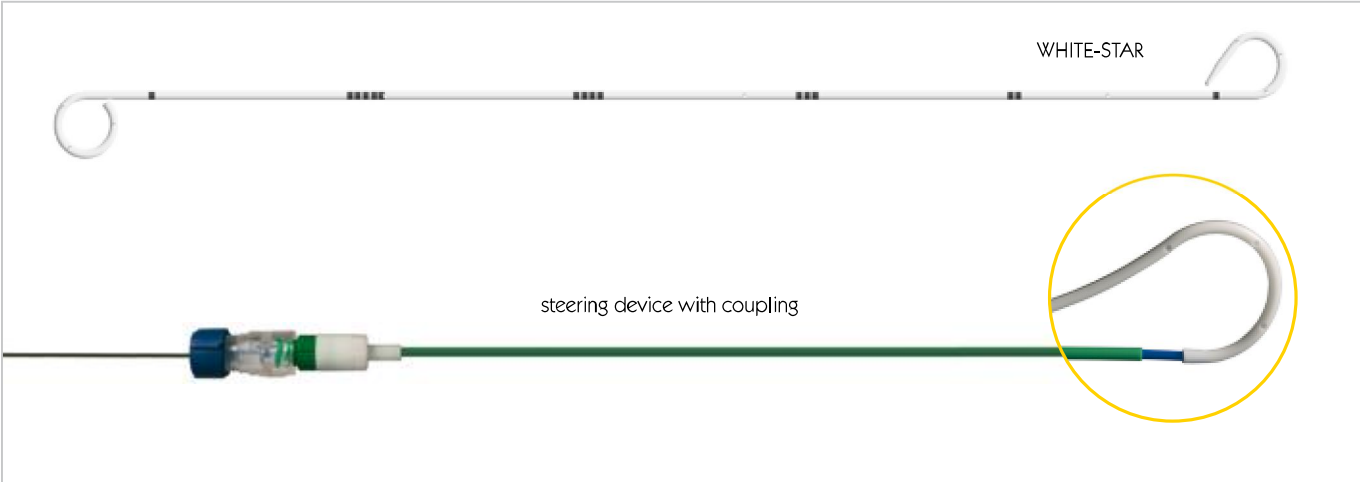

Ureteral Stent



- ◆ Polyurethane (PUR)
- ◆ Drainage openings in pigtails and shaft
- ◆ Open or closed tip
- ◆ Markings for exact positioning
- ◆ Radiopaque

Steering Device with Coupling

- ◆ Pusher
- ◆ Disconnecter
- ◆ Coupling to fasten guide wire
- ◆ Length: 43 cm

Guide Wire

- ◆ PTFE-coated

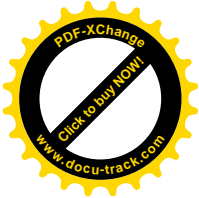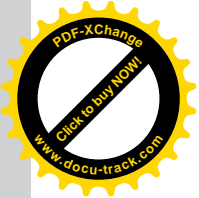

WHITE-STAR Double-J Ureteral Stent-Set


Steerable with coupling

Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24	ES-570524	ES-580524	145	0.028"	straight	standard
4.7	26	ES-570526	ES-580526	145	0.028"	straight	standard
06	24	ES-570624	ES-580624	150	0.035"	straight	standard / movable core
06	26	ES-570626	ES-580626	150	0.035"	straight	standard / movable core
06	28	ES-570628	ES-580628	150	0.035"	straight	standard / movable core
06	30	ES-570630	ES-580630	150	0.035"	straight	standard / movable core
07	24	ES-570724	ES-580724	150	0.035"	straight	standard / movable core
07	26	ES-570726	ES-580726	150	0.035"	straight	standard / movable core
07	28	ES-570728	ES-580728	150	0.035"	straight	standard / movable core
07	30	ES-570730	ES-580730	150	0.035"	straight	standard / movable core
08	24	ES-570824	ES-580824	150	0.035"	straight	standard / movable core
08	26	ES-570826	ES-580826	150	0.035"	straight	standard / movable core
08	28	ES-570828	ES-580828	150	0.035"	straight	standard / movable core
08	30	ES-570830	ES-580830	150	0.035"	straight	standard / movable core

Closed tip							
Size (FR/CH)	Length (cm)		Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24		ES-590524	125	0.028"	straight	standard
4.7	26		ES-590526	125	0.028"	straight	standard
06	24		ES-590624	125	0.035"	J-tip	standard / movable core
06	26		ES-590626	125	0.035"	J-tip	standard / movable core
06	28		ES-590628	125	0.035"	J-tip	standard / movable core
06	30		ES-590630	125	0.035"	J-tip	standard / movable core
07	24		ES-590724	125	0.035"	J-tip	standard / movable core
07	26		ES-590726	125	0.035"	J-tip	standard / movable core
07	28		ES-590728	125	0.035"	J-tip	standard / movable core
07	30		ES-590730	125	0.035"	J-tip	standard / movable core
08	24		ES-590824	125	0.035"	J-tip	standard / movable core
08	26		ES-590826	125	0.035"	J-tip	standard / movable core
08	28		ES-590828	125	0.035"	J-tip	standard / movable core
08	30		ES-590830	125	0.035"	J-tip	standard / movable core

WHITE-STAR MONO Mono-J Ureteral Stent-Set

SET COMPONENTS

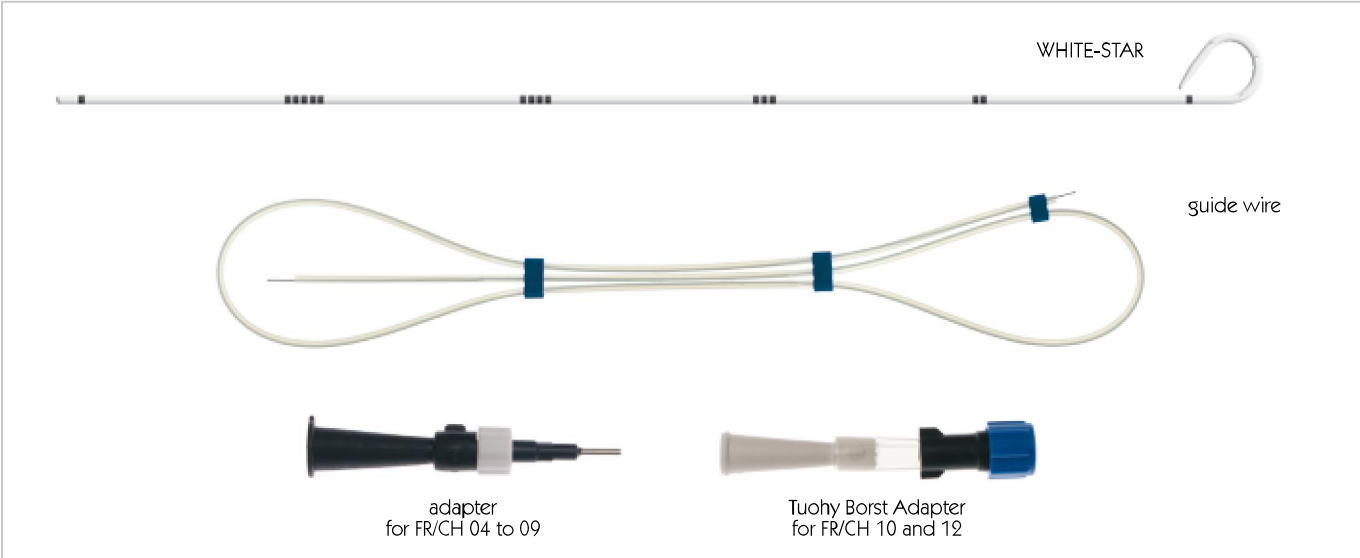

Ureteral Stent

- ♦ Polyurethane (PUR)
- ♦ Drainage openings only in pigtail
- ♦ Open or closed tip
- ♦ Markings for exact positioning
- ♦ Radiopaque

Guide Wire

- ♦ PTFE-coated



Adapter

- ♦ For attachment of the urinary drainage bag

Clamp

- ♦ Only in sets with closed tip



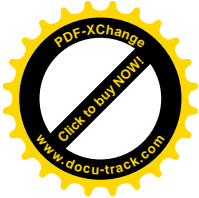
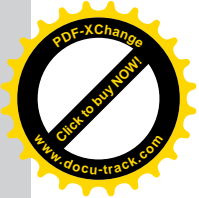


URETERAL STENTS

WHITE-STAR MONO Mono-J Ureteral Stent-Set

Open tip						
Size (FR/CH)	Length (cm)	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
04	70	EN-350471	175	0.028"	straight	standard
4.7	70	EN-350571	175	0.028"	straight	standard
06	70	EN-350671	175	0.035"	straight	standard
07	70	EN-350771	175	0.035"	straight	standard
08	70	EN-350871	175	0.035"	straight	standard
09	70	EN-350971	175	0.045"	straight	standard
10	70	EN-351071	175	0.045"	straight	standard
12	70	EN-351271	175	0.045"	straight	standard


Closed tip						
Size (FR/CH)	Length (cm)	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
04	70	EN-350470	125	0.028"	straight	standard
4.7	70	EN-350570	125	0.028"	straight	standard
06	70	EN-350670	125	0.035"	straight	standard
07	70	EN-350770	125	0.035"	straight	standard
08	70	EN-350870	125	0.035"	straight	standard
09	70	EN-350970	125	0.045"	straight	standard
10	70	EN-351070	125	0.045"	straight	standard
12	70	EN-351270	125	0.045"	straight	standard

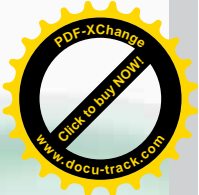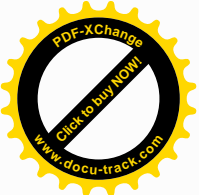

URETERAL STENTS

GREEN-STAR




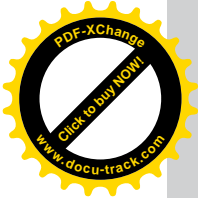

GREEN-STAR Double-J Ureteral Stents-Set

Soft material

SET COMPONENTS

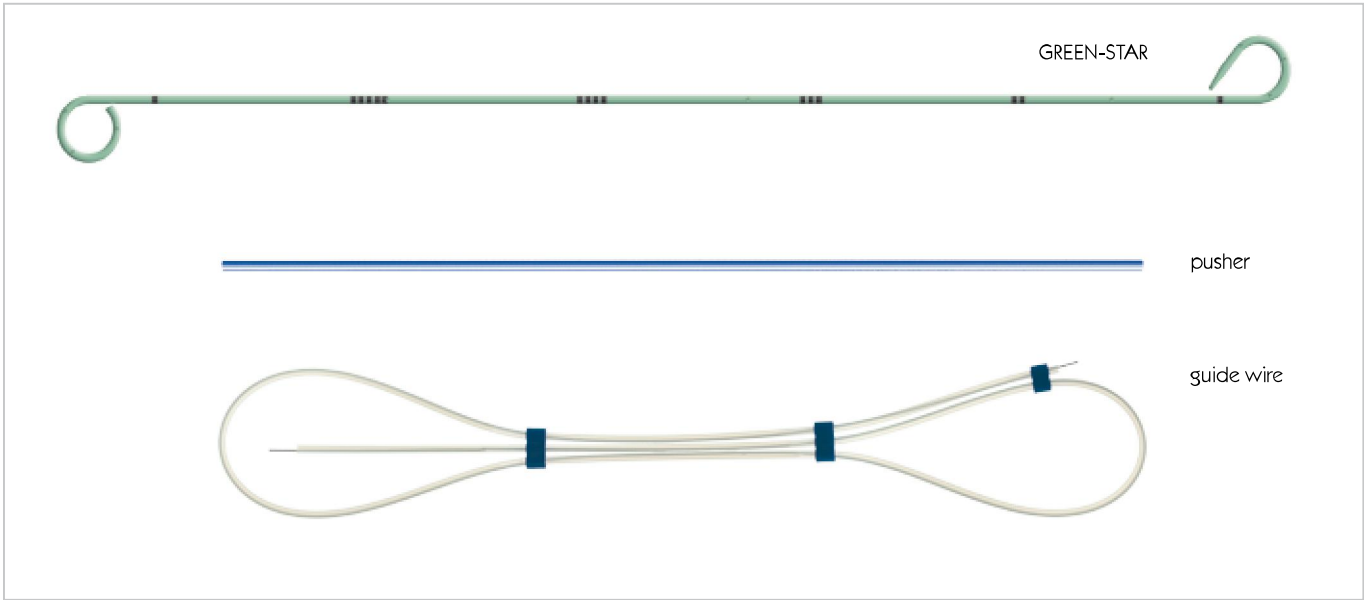

Ureteral Stent

- ◆ Polyurethane (PUR)
- ◆ Special soft material for high patient comfort
- ◆ Drainage openings in pigtails and shaft
- ◆ Open or closed tip
- ◆ Markings for exact positioning
- ◆ High radiopaque visibility

Pusher

- ◆ Length: 40 cm



Guide Wire

- ◆ PTFE-coated

Clamp

- ◆ Only in sets with closed tip



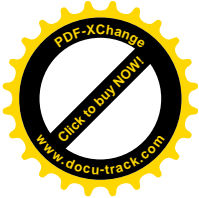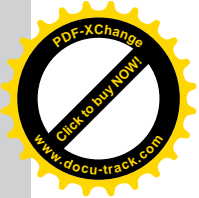


URETERAL STENTS

GREEN-STAR Double-J Ureteral Stents-Set

Soft material

Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	24	EG-410624	EG-360624	125	0.035"	straight	standard
06	26	EG-410626	EG-360626	125	0.035"	straight	standard
06	28	EG-410628	EG-360628	125	0.035"	straight	standard
06	30	EG-410630	EG-360630	125	0.035"	straight	standard
07	24	EG-410724	EG-360724	125	0.035"	straight	standard
07	26	EG-410726	EG-360726	125	0.035"	straight	standard
07	28	EG-410728	EG-360728	125	0.035"	straight	standard
07	30	EG-410730	EG-360730	125	0.035"	straight	standard
08	24	EG-410824	EG-360824	125	0.035"	straight	standard
08	26	EG-410826	EG-360826	125	0.035"	straight	standard
08	28	EG-410828	EG-360828	125	0.035"	straight	standard
08	30	EG-410830	EG-360830	125	0.035"	straight	standard


Closed tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	24	EG-400624	EG-390624	125	0.035"	straight	standard
06	26	EG-400626	EG-390626	125	0.035"	straight	standard
06	28	EG-400628	EG-390628	125	0.035"	straight	standard
06	30	EG-400630	EG-390630	125	0.035"	straight	standard
07	24	EG-400724	EG-390724	125	0.035"	straight	standard
07	26	EG-400726	EG-390726	125	0.035"	straight	standard
07	28	EG-400728	EG-390728	125	0.035"	straight	standard
07	30	EG-400730	EG-390730	125	0.035"	straight	standard
08	24	EG-400824	EG-390824	125	0.035"	straight	standard
08	26	EG-400826	EG-390826	125	0.035"	straight	standard
08	28	EG-400828	EG-390828	125	0.035"	straight	standard
08	30	EG-400830	EG-390830	125	0.035"	straight	standard

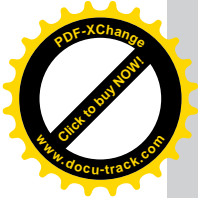

GREEN-STAR Double-J Ureteral Stents-Set


Soft material, flexible length

SET COMPONENTS

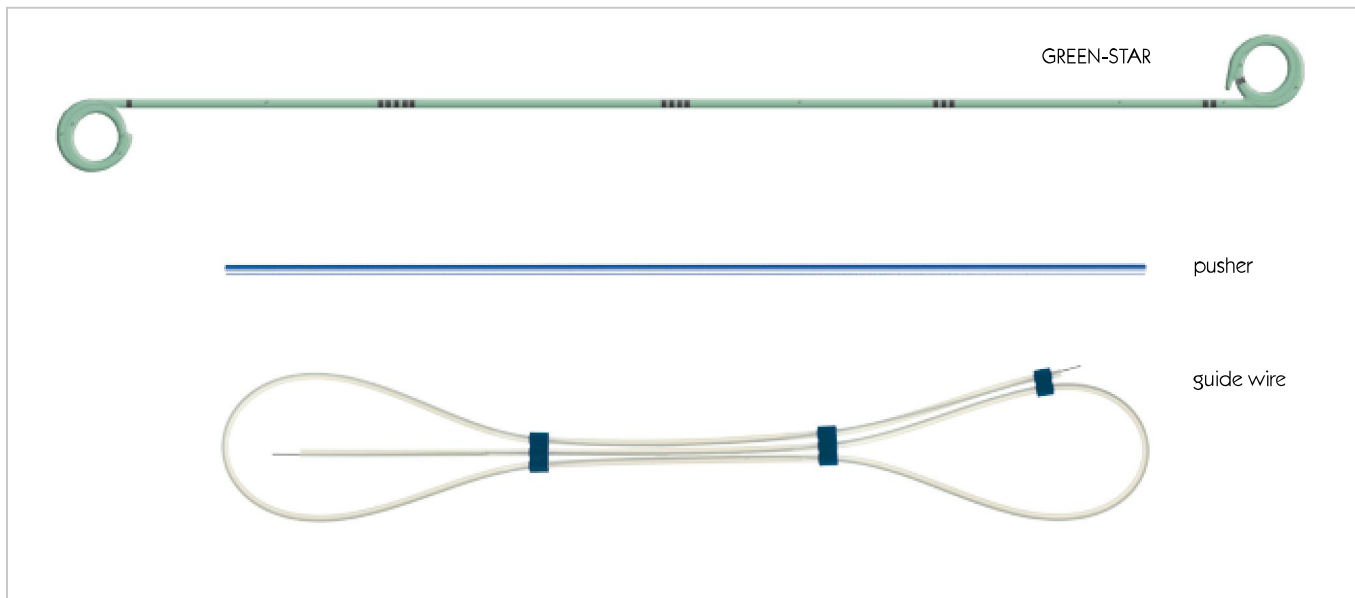

Ureteral Stent

- ◆ Polyurethane (PUR)
- ◆ Special soft material for high patient comfort
- ◆ Multiple coils on both ends for flexible shaft length from 22 to 30 cm
- ◆ Drainage openings in pigtails and shaft
- ◆ Open or closed tip
- ◆ Markings for exact positioning
- ◆ High radiopaque visibility

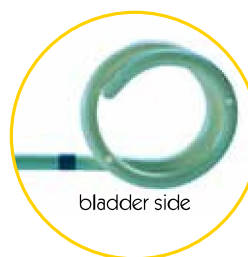

Pusher

- ◆ Length: 40 cm



Guide Wire

- ◆ PTFE-coated

Clamp

- ◆ Only in sets with closed tip



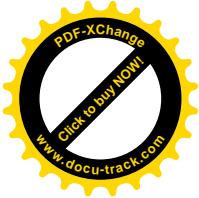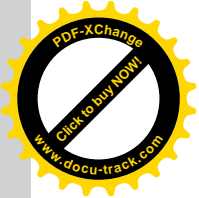



URETERAL STENTS

GREEN-STAR Double-J Ureteral Stents-Set

Soft material, flexible length

Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	22-30	EG-410622-30-BS	EG-360622-30-BS	125	0.035"	straight	standard
07	22-30	EG-410722-30-BS	EG-360722-30-BS	125	0.035"	straight	standard
08	22-30	EG-410822-30-BS	EG-360822-30-BS	125	0.035"	straight	standard

Closed tip							
Size (FR/CH)	Length (cm)		Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	22-30		EG-390622-30-BS	125	0.035"	straight	standard
07	22-30		EG-390722-30-BS	125	0.035"	straight	standard
08	22-30		EG-390822-30-BS	125	0.035"	straight	standard

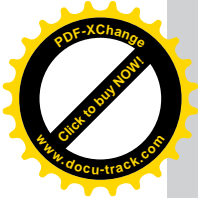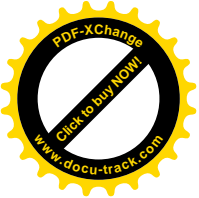

GREEN-STAR Double-J Ureteral Stents-Set

Soft material, steerable

SET COMPONENTS



Ureteral Stent

- ◆ Polyurethane (PUR)
- ◆ Special soft material for high patient comfort
- ◆ Drainage openings in pigtails and shaft
- ◆ Open or closed tip
- ◆ Markings for exact positioning
- ◆ High radiopaque visibility

Steering Device

- ◆ Pusher
- ◆ Disconnecter
- ◆ Length: 53 cm



Guide Wire

- ◆ PTFE-coated

Clamp

- ◆ Only in sets with closed tip



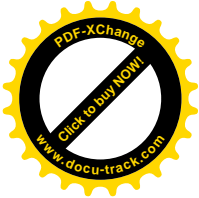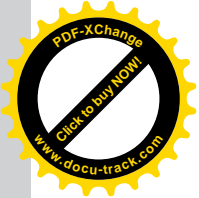

URETERAL STENTS

GREEN-STAR Double-J Ureteral Stents-Set

Soft material, steerable


Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	24	EG-470624	EG-480624	150	0.035"	straight	standard / movable core
06	26	EG-470626	EG-480626	150	0.035"	straight	standard / movable core
06	28	EG-470628	EG-480628	150	0.035"	straight	standard / movable core
06	30	EG-470630	EG-480630	150	0.035"	straight	standard / movable core
07	24	EG-470724	EG-480724	150	0.035"	straight	standard / movable core
07	26	EG-470726	EG-480726	150	0.035"	straight	standard / movable core
07	28	EG-470728	EG-480728	150	0.035"	straight	standard / movable core
07	30	EG-470730	EG-480730	150	0.035"	straight	standard / movable core
08	24	EG-470824	EG-480824	150	0.035"	straight	standard / movable core
08	26	EG-470826	EG-480826	150	0.035"	straight	standard / movable core
08	28	EG-470828	EG-480828	150	0.035"	straight	standard / movable core
08	30	EG-470830	EG-480830	150	0.035"	straight	standard / movable core

Closed tip							
Size (FR/CH)	Length (cm)		Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	24		EG-490624	125	0.035"	J-tip	standard / movable core
06	26		EG-490626	125	0.035"	J-tip	standard / movable core
06	28		EG-490628	125	0.035"	J-tip	standard / movable core
06	30		EG-490630	125	0.035"	J-tip	standard / movable core
07	24		EG-490724	125	0.035"	J-tip	standard / movable core
07	26		EG-490726	125	0.035"	J-tip	standard / movable core
07	28		EG-490728	125	0.035"	J-tip	standard / movable core
07	30		EG-490730	125	0.035"	J-tip	standard / movable core
08	24		EG-490824	125	0.035"	J-tip	standard / movable core
08	26		EG-490826	125	0.035"	J-tip	standard / movable core
08	28		EG-490828	125	0.035"	J-tip	standard / movable core
08	30		EG-490830	125	0.035"	J-tip	standard / movable core

GREEN-STAR Double-J Ureteral Stents-Set

Soft material, steerable with coupling

SET COMPONENTS

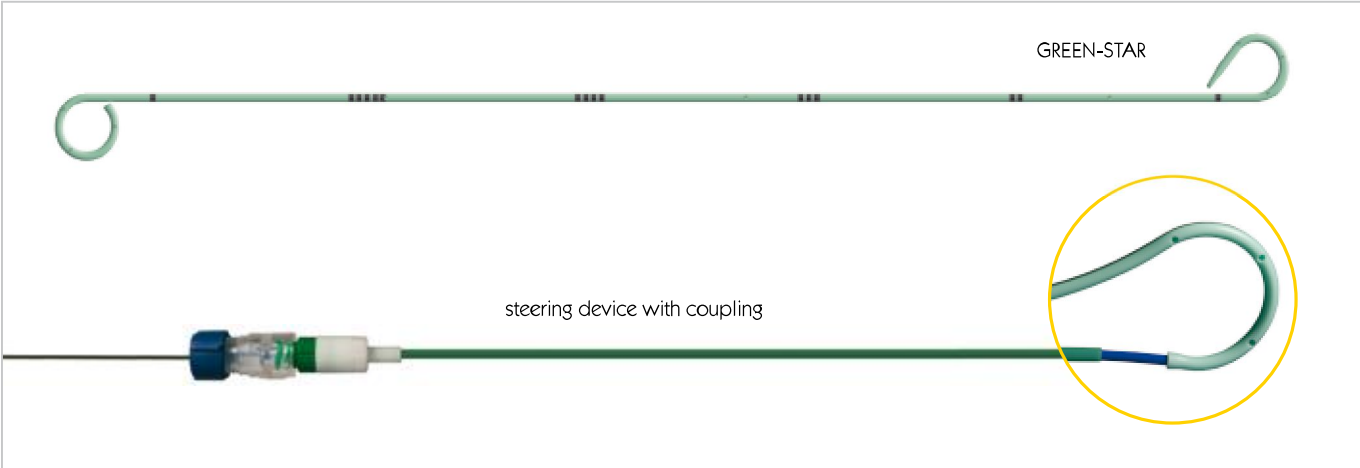

Ureteral Stent


- ◆ Polyurethane (PUR)
- ◆ Special soft material for high patient comfort
- ◆ Drainage openings in pigtails and shaft
- ◆ Open or closed tip
- ◆ Markings for exact positioning
- ◆ High radiopaque visibility

Steering Device with Coupling

- ◆ Pusher
- ◆ Disconnecter
- ◆ Coupling to fasten guide wire
- ◆ Length: 43 cm

Guide Wire

- ◆ PTFE-coated

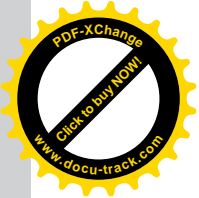

GREEN-STAR Double-J Ureteral Stents-Set

Soft material, steerable with coupling


Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	24	EG-570624	EG-580624	150	0.035"	straight	standard / movable core
06	26	EG-570626	EG-580626	150	0.035"	straight	standard / movable core
06	28	EG-570628	EG-580628	150	0.035"	straight	standard / movable core
06	30	EG-570630	EG-580630	150	0.035"	straight	standard / movable core
07	24	EG-570724	EG-580724	150	0.035"	straight	standard / movable core
07	26	EG-570726	EG-580726	150	0.035"	straight	standard / movable core
07	28	EG-570728	EG-580728	150	0.035"	straight	standard / movable core
07	30	EG-570730	EG-580730	150	0.035"	straight	standard / movable core
08	24	EG-570824	EG-580824	150	0.035"	straight	standard / movable core
08	26	EG-570826	EG-580826	150	0.035"	straight	standard / movable core
08	28	EG-570828	EG-580828	150	0.035"	straight	standard / movable core
08	30	EG-570830	EG-580830	150	0.035"	straight	standard / movable core

Closed tip							
Size (FR/CH)	Length (cm)		Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	24		EG-590624	125	0.035"	J-tip	standard / movable core
06	26		EG-590626	125	0.035"	J-tip	standard / movable core
06	28		EG-590628	125	0.035"	J-tip	standard / movable core
06	30		EG-590630	125	0.035"	J-tip	standard / movable core
07	24		EG-590724	125	0.035"	J-tip	standard / movable core
07	26		EG-590726	125	0.035"	J-tip	standard / movable core
07	28		EG-590728	125	0.035"	J-tip	standard / movable core
07	30		EG-590730	125	0.035"	J-tip	standard / movable core
08	24		EG-590824	125	0.035"	J-tip	standard / movable core
08	26		EG-590826	125	0.035"	J-tip	standard / movable core
08	28		EG-590828	125	0.035"	J-tip	standard / movable core
08	30		EG-590830	125	0.035"	J-tip	standard / movable core

GREEN-STAR MONO Mono-J Ureteral Stents-Set

Soft material

SET COMPONENTS

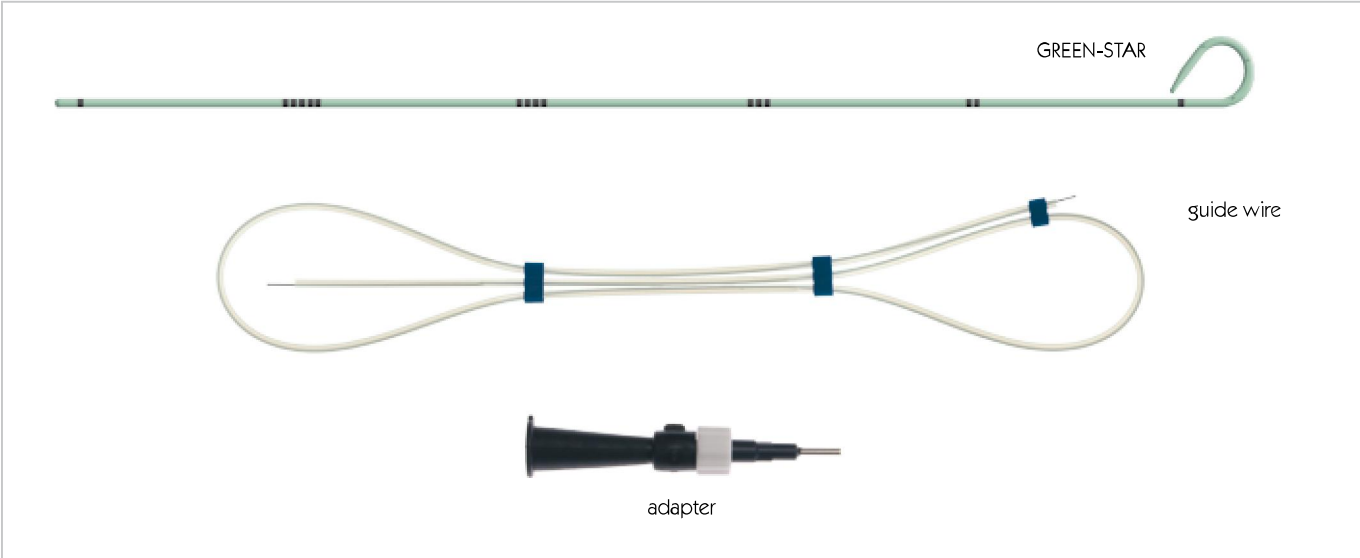

Ureteral Stent

- ♦ Polyurethane (PUR)
- ♦ Special soft material for high patient comfort
- ♦ Drainage openings only in pigtail
- ♦ Open or closed tip
- ♦ Markings for exact positioning
- ♦ High radiopaque visibility

Guide Wire

- ♦ PTFE-coated



Adapter

- ♦ For attachment of the urinary drainage bag

Clamp

- ♦ Only in sets with closed tip



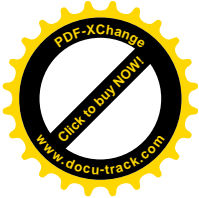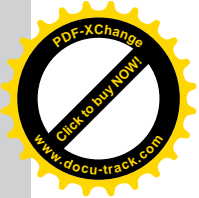



URETERAL STENTS

GREEN-STAR MONO Mono-J Ureteral Stents-Set

Soft material

Open tip						
Size (FR/CH)	Length (cm)	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	70	EG-350671	175	0.035"	straight	standard
07	70	EG-350771	175	0.035"	straight	standard
08	70	EG-350871	175	0.035"	straight	standard

Closed tip						
Size (FR/CH)	Length (cm)	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	70	EG-350670	125	0.035"	straight	standard
07	70	EG-350770	125	0.035"	straight	standard
08	70	EG-350870	125	0.035"	straight	standard

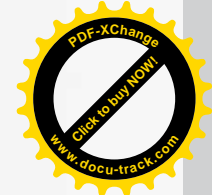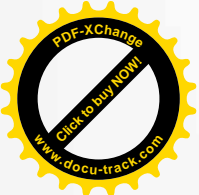

URETERAL STENTS

SPECIAL STENTS

YELLOW-STAR TUMOR-STENT

WHITE-STAR ESWL

WHITE-STAR INTRAOPERATIVE

WHITE-STAR STENOSIS


URETERAL STENT for supravescical urine drainage



YELLOW-STAR TUMOR-STENT Double-J Ureteral Stents-Set

PC-coated

SET COMPONENTS

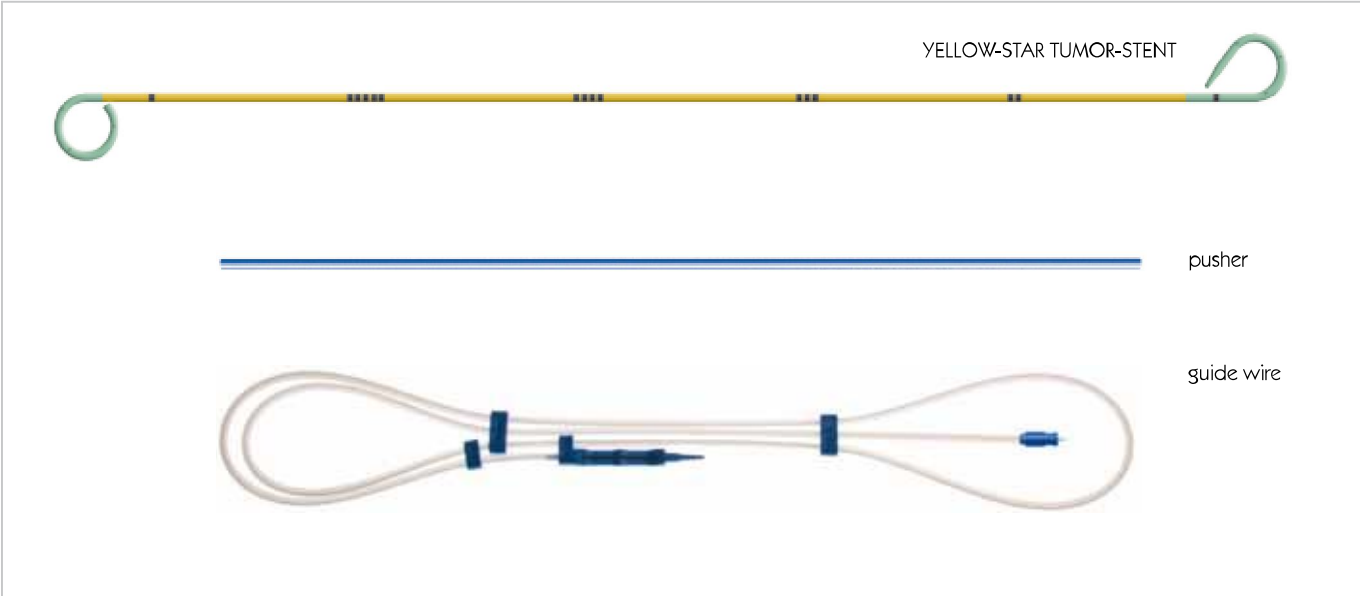

Ureteral Stent



- ♦ Polyurethane (PUR)
- ♦ PC coating to minimise risk of infection and encrustation, resulting in extended indwelling time
- ♦ Pressure resistant shaft, without drainage openings for preventing ingrowth of tumor tissue
- ♦ Soft pigtails for better wearing comfort
- ♦ Drainage openings only in pigtails
- ♦ Open tip
- ♦ Markings for exact positioning
- ♦ Radiopaque

Pusher

- ♦ Length: 40 cm

Guide Wire

- ♦ Hydrophilic coating for better gliding ability – prior to application, guide wire must be moistened in dispenser

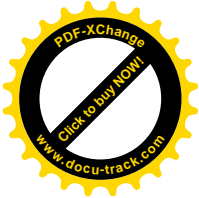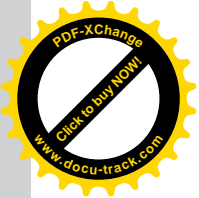

URETERAL STENTS


YELLOW-STAR TUMOR-STENT Double-J Ureteral Stents-Set PC-coated

Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	22	TU-410622	TU-360622	125	0.035"	straight	standard
06	24	TU-410624	TU-360624	125	0.035"	straight	standard
06	26	TU-410626	TU-360626	125	0.035"	straight	standard
06	28	TU-410628	TU-360628	125	0.035"	straight	standard
06	30	TU-410630	TU-360630	125	0.035"	straight	standard
06	32	TU-410632	TU-360632	125	0.035"	straight	standard
07	22	TU-410722	TU-360722	125	0.035"	straight	standard
07	24	TU-410724	TU-360724	125	0.035"	straight	standard
07	26	TU-410726	TU-360726	125	0.035"	straight	standard
07	28	TU-410728	TU-360728	125	0.035"	straight	standard
07	30	TU-410730	TU-360730	125	0.035"	straight	standard
07	32	TU-410732	TU-360732	125	0.035"	straight	standard
08	22	TU-410822	TU-360822	125	0.035"	straight	standard
08	24	TU-410824	TU-360824	125	0.035"	straight	standard
08	26	TU-410826	TU-360826	125	0.035"	straight	standard
08	28	TU-410828	TU-360828	125	0.035"	straight	standard
08	30	TU-410830	TU-360830	125	0.035"	straight	standard
08	32	TU-410832	TU-360832	125	0.035"	straight	standard

YELLOW-STAR TUMOR-STENT Double-J Ureteral Stents-Set

PC-coated, steerable

SET COMPONENTS

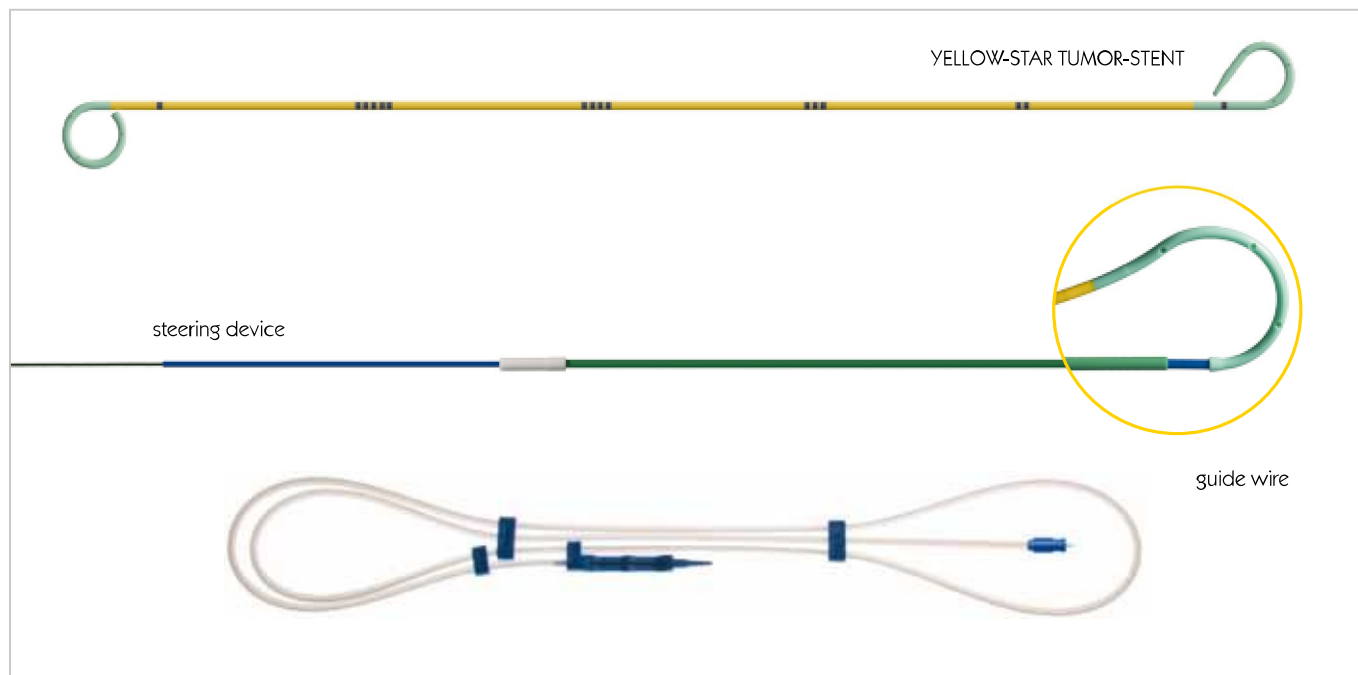

Ureteral Stent



- ◆ Polyurethane (PUR)
- ◆ PC coating to minimise risk of infection and encrustation, resulting in extended indwelling time
- ◆ Pressure resistant shaft, without drainage openings for preventing ingrowth of tumor tissue
- ◆ Soft pigtails for better wearing comfort
- ◆ Drainage openings only in pigtails
- ◆ Open tip
- ◆ Markings for exact positioning
- ◆ Radiopaque

Steering Device

- ◆ Pusher
- ◆ Disconnecter
- ◆ Length: 53 cm

Guide Wire

- ◆ Hydrophilic coating for better gliding ability – prior to application, guide wire must be moistened in dispenser

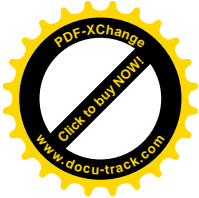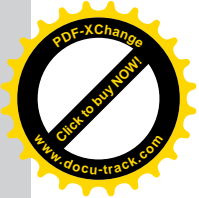

URETERAL STENTS

YELLOW-STAR TUMOR-STENT Double-J Ureteral Stents-Set

PC-coated, steerable

Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	22	TU-470622	TU-480622	150	0.035"	straight	standard
06	24	TU-470624	TU-480624	150	0.035"	straight	standard
06	26	TU-470626	TU-480626	150	0.035"	straight	standard
06	28	TU-470628	TU-480628	150	0.035"	straight	standard
06	30	TU-470630	TU-480630	150	0.035"	straight	standard
06	32	TU-470632	TU-480632	150	0.035"	straight	standard
07	22	TU-470722	TU-480722	150	0.035"	straight	standard
07	24	TU-470724	TU-480724	150	0.035"	straight	standard
07	26	TU-470726	TU-480726	150	0.035"	straight	standard
07	28	TU-470728	TU-480728	150	0.035"	straight	standard
07	30	TU-470730	TU-480730	150	0.035"	straight	standard
07	32	TU-470732	TU-480732	150	0.035"	straight	standard
08	22	TU-470822	TU-480822	150	0.035"	straight	standard
08	24	TU-470824	TU-480824	150	0.035"	straight	standard
08	26	TU-470826	TU-480826	150	0.035"	straight	standard
08	28	TU-470828	TU-480828	150	0.035"	straight	standard
08	30	TU-470830	TU-480830	150	0.035"	straight	standard
08	32	TU-470832	TU-480832	150	0.035"	straight	standard

YELLOW-STAR TUMOR-STENT Double-J Ureteral Stent-Set PC-coated, steerable with coupling

SET COMPONENTS



Ureteral Stent


- ◆ Polyurethane (PUR)
- ◆ PC coating to minimise risk of infection and encrustation, resulting in extended indwelling time
- ◆ Pressure resistant shaft, without drainage openings for preventing ingrowth of tumor tissue
- ◆ Soft pigtails for better wearing comfort
- ◆ Drainage openings only in pigtails
- ◆ Open tip
- ◆ Markings for exact positioning
- ◆ Radiopaque

Steering Device with Coupling

- ◆ Pusher
- ◆ Disconnecter
- ◆ Coupling to fasten guide wire
- ◆ Length: 43 cm

Guide Wire

- ◆ Hydrophilic coating for better gliding ability – prior to application, guide wire must be moistened in dispenser

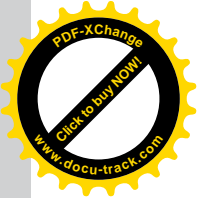


URETERAL STENTS

YELLOW-STAR TUMOR-STENT Double-J Ureteral Stent-Set

PC-coated, steerable with coupling


Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
06	22	TU-570622	TU-580622	150	0.035"	straight	standard
06	24	TU-570624	TU-580624	150	0.035"	straight	standard
06	26	TU-570626	TU-580626	150	0.035"	straight	standard
06	28	TU-570628	TU-580628	150	0.035"	straight	standard
06	30	TU-570630	TU-580630	150	0.035"	straight	standard
06	32	TU-570632	TU-580632	150	0.035"	straight	standard
07	22	TU-570722	TU-580722	150	0.035"	straight	standard
07	24	TU-570724	TU-580724	150	0.035"	straight	standard
07	26	TU-570726	TU-580726	150	0.035"	straight	standard
07	28	TU-570728	TU-580728	150	0.035"	straight	standard
07	30	TU-570730	TU-580730	150	0.035"	straight	standard
07	32	TU-570732	TU-580732	150	0.035"	straight	standard
08	22	TU-570822	TU-580822	150	0.035"	straight	standard
08	24	TU-570824	TU-580824	150	0.035"	straight	standard
08	26	TU-570826	TU-580826	150	0.035"	straight	standard
08	28	TU-570828	TU-580828	150	0.035"	straight	standard
08	30	TU-570830	TU-580830	150	0.035"	straight	standard
08	32	TU-570832	TU-580832	150	0.035"	straight	standard

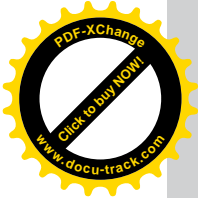

WHITE-STAR ESWL Double-J Ureteral Stents-Set


With special profile

SET COMPONENTS

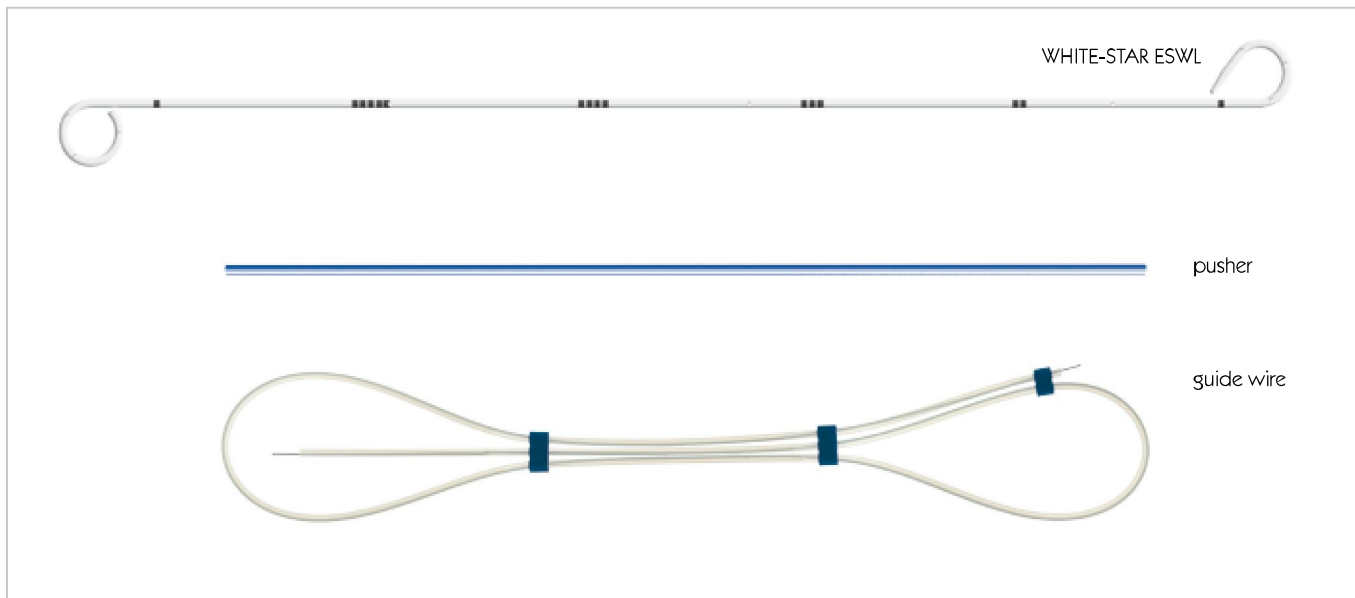

Ureteral Stent

- ♦ Polyurethane (PUR)
- ♦ With special profile for supporting transport of stone debris
- ♦ Drainage openings only in pigtails
- ♦ Open or closed tip
- ♦ Markings for exact positioning
- ♦ Radiopaque

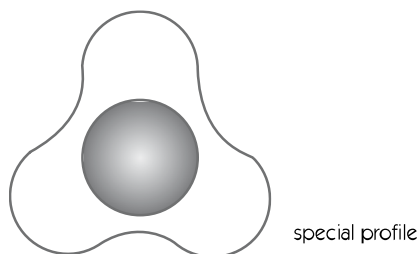

Pusher

- ♦ Length: 40 cm



Guide Wire

- ♦ PTFE-coated

Clamp

- ♦ Only in sets with closed tip



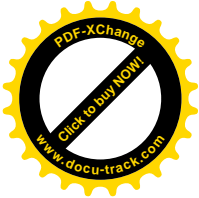
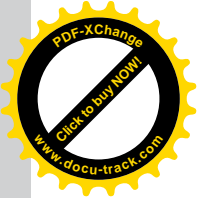


URETERAL STENTS

WHITE-STAR ESWL Double-J Ureteral Stents-Set With special profile

Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
07	24	EL-410724	EL-360724	125	0.035"	straight	standard
07	26	EL-410726	EL-360726	125	0.035"	straight	standard
07	28	EL-410728	EL-360728	125	0.035"	straight	standard
07	30	EL-410730	EL-360730	125	0.035"	straight	standard


Closed tip							
Size (FR/CH)	Length (cm)		Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
07	24		EL-390724	125	0.035"	straight	standard
07	26		EL-390726	125	0.035"	straight	standard
07	28		EL-390728	125	0.035"	straight	standard
07	30		EL-390730	125	0.035"	straight	standard

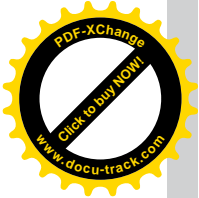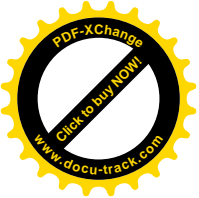

WHITE-STAR ESWL Double-J Ureteral Stent-Set


With special profile, steerable

SET COMPONENTS

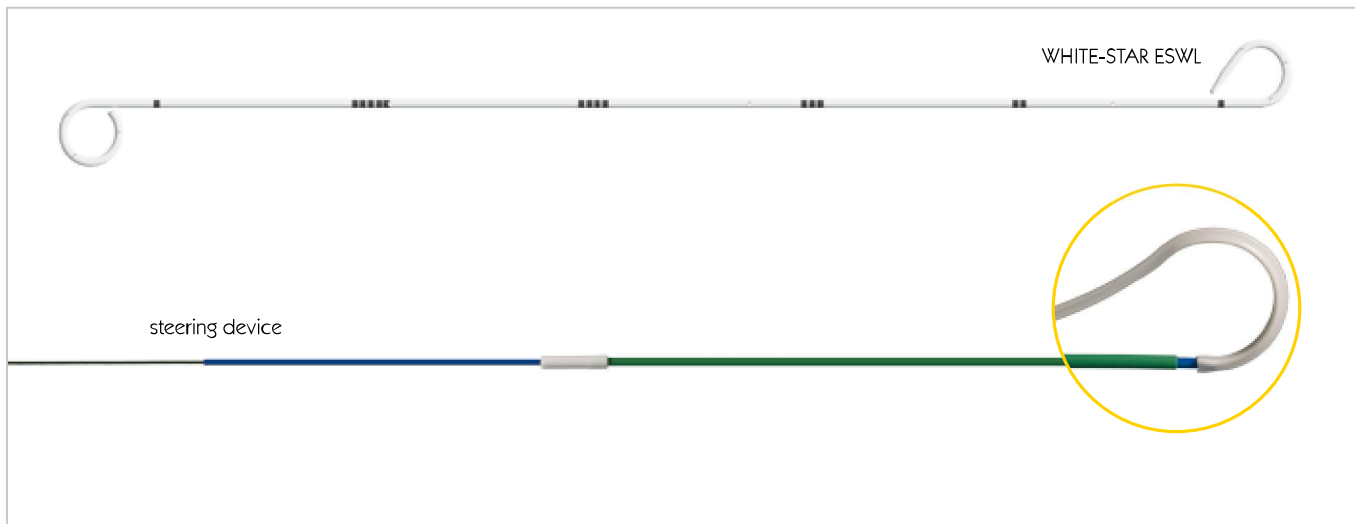

Ureteral Stent

- ♦ Polyurethane (PUR)
- ♦ With special profile for supporting transport of stone debris
- ♦ Drainage openings only in pigtails
- ♦ Open or closed tip
- ♦ Markings for exact positioning
- ♦ Radiopaque

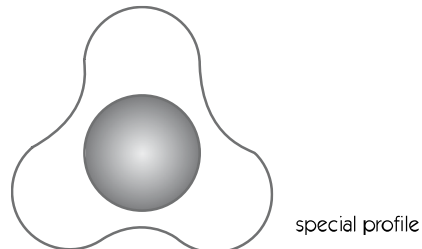

Steering Device

- ♦ Pusher
- ♦ Disconnecter
- ♦ Length: 53 cm



Guide Wire

- ♦ PTFE-coated

Clamp

- ♦ Only in sets with closed tip



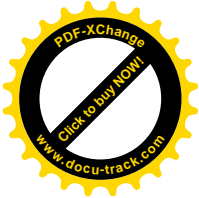
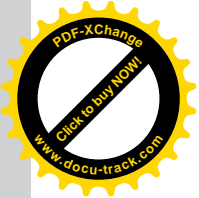


URETERAL STENTS

WHITE-STAR ESWL Double-J Ureteral Stent-Set With special profile, steerable

Open tip							
Size (FR/CH)	Length (cm)	Order Code without Guide Wire	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
07	24	EL-470724	EL-480724	150	0.035"	straight	standard / movable core
07	26	EL-470726	EL-480726	150	0.035"	straight	standard / movable core
07	28	EL-470728	EL-480728	150	0.035"	straight	standard / movable core
07	30	EL-470730	EL-480730	150	0.035"	straight	standard / movable core


Closed tip							
Size (FR/CH)	Length (cm)		Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
07	24		EL-490724	125	0.035"	J-Tip	standard / movable core
07	26		EL-490726	125	0.035"	J-Tip	standard / movable core
07	28		EL-490728	125	0.035"	J-Tip	standard / movable core
07	30		EL-490730	125	0.035"	J-Tip	standard / movable core

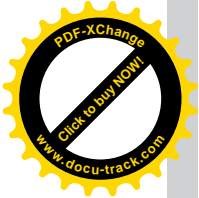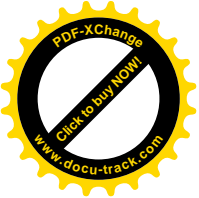

WHITE-STAR INTRAOPERATIVE Double-J Ureteral Stents-Set

Round tip

SET COMPONENTS

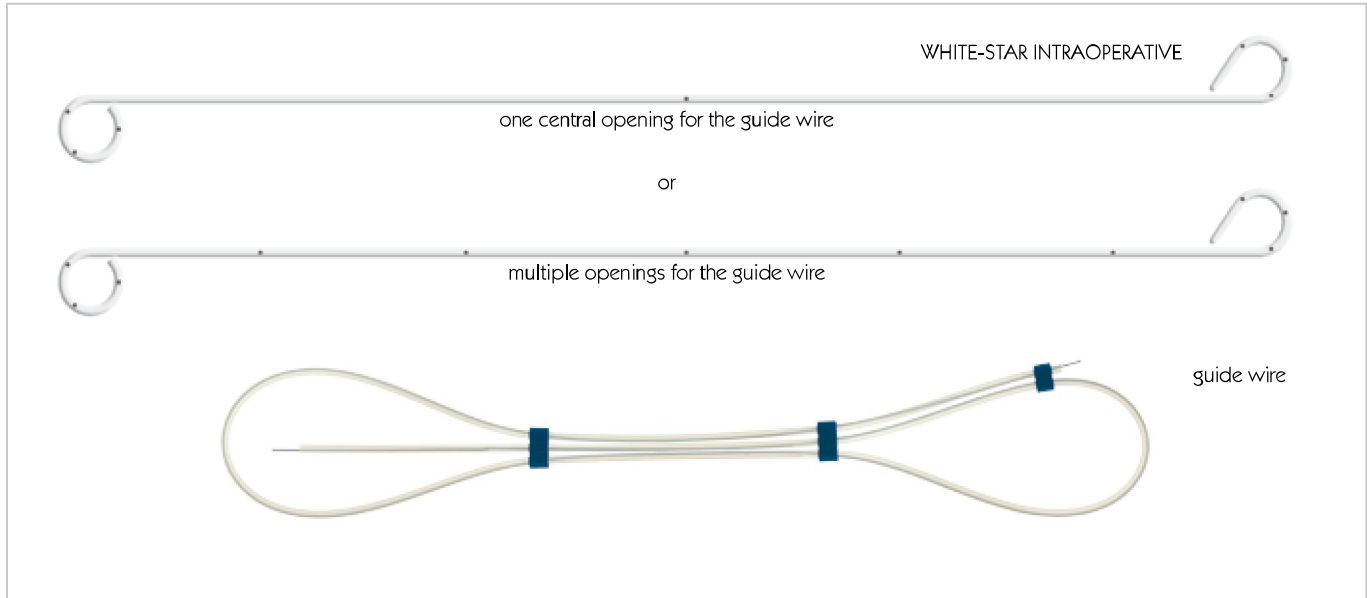

Ureteral Stent

- ◆ Polyurethane (PUR)
- ◆ Both ends closed
- ◆ At the shaft one central opening or multiple openings for the guide wire
- ◆ Radiopaque





bladder side



kidney side

Guide Wire

- ◆ PTFE-coated

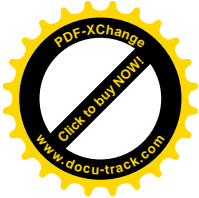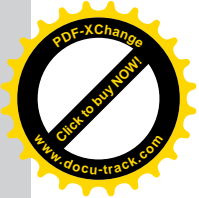

WHITE-STAR INTRAOPERATIVE Double-J Ureteral Stents-Set


Round tip

One central opening for the guide wire – closed tip

Size (FR/CH)	Length (cm)	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24	EN-320524	100	0.028"	straight	standard
4.7	26	EN-320526	100	0.028"	straight	standard
4.7	28	EN-320528	100	0.028"	straight	standard
4.7	30	EN-320530	100	0.028"	straight	standard
06	24	EN-320624	100	0.035"	straight	standard
06	26	EN-320626	100	0.035"	straight	standard
06	28	EN-320628	100	0.035"	straight	standard
06	30	EN-320630	100	0.035"	straight	standard
07	24	EN-320724	100	0.035"	straight	standard
07	26	EN-320726	100	0.035"	straight	standard
07	28	EN-320728	100	0.035"	straight	standard
07	30	EN-320730	100	0.035"	straight	standard
08	24	EN-320824	100	0.035"	straight	standard
08	26	EN-320826	100	0.035"	straight	standard
08	28	EN-320828	100	0.035"	straight	standard
08	30	EN-320830	100	0.035"	straight	standard

Multiple openings for the guide wire – closed tip

Size (FR/CH)	Length (cm)	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
4.7	24	EN-330524	100	0.028"	straight	standard
4.7	26	EN-330526	100	0.028"	straight	standard
4.7	28	EN-330528	100	0.028"	straight	standard
4.7	30	EN-330530	100	0.028"	straight	standard
06	24	EN-330624	100	0.035"	straight	standard
06	26	EN-330626	100	0.035"	straight	standard
06	28	EN-330628	100	0.035"	straight	standard
06	30	EN-330630	100	0.035"	straight	standard
07	24	EN-330724	100	0.035"	straight	standard
07	26	EN-330726	100	0.035"	straight	standard
07	28	EN-330728	100	0.035"	straight	standard
07	30	EN-330730	100	0.035"	straight	standard
08	24	EN-330824	100	0.035"	straight	standard
08	26	EN-330826	100	0.035"	straight	standard
08	28	EN-330828	100	0.035"	straight	standard
08	30	EN-330830	100	0.035"	straight	standard

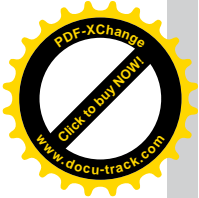

WHITE-STAR STENOSIS Stenosis Stent-Set


Partial shaft enlargement

SET COMPONENTS



Ureteral Stent

- ♦ Polyurethane (PUR)
- ♦ Drainage openings only in pigtails
- ♦ Open tip
- ♦ Markings for exact positioning
- ♦ Radiopaque
- ♦ Partial enlargement (proximal, central, distal) for dilation of stenosis

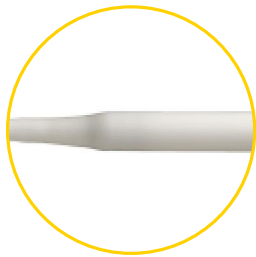


Pusher

- ♦ Length: 40 cm

Guide Wire

- ♦ PTFE-coated

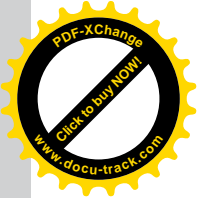

WHITE-STAR STENOSIS Stenosis Stent-Set

Partial shaft enlargement

Bladder side / proximal enlargement – open tip


Size (FR/CH)	Shaft Enlargement (FR/CH)	Length (cm)	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
07	12	26	TU-201226-P	125	0.035"	straight	standard
07	12	28	TU-201228-P	125	0.035"	straight	standard
07	12	30	TU-201230-P	125	0.035"	straight	standard

Central enlargement – open tip

Size (FR/CH)	Shaft Enlargement (FR/CH)	Length (cm)	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
07	12	26	TU-201226-M	125	0.035"	straight	standard
07	12	28	TU-201228-M	125	0.035"	straight	standard
07	12	30	TU-201230-M	125	0.035"	straight	standard

Kidney side / distal enlargement – open tip

Size (FR/CH)	Shaft Enlargement (FR/CH)	Length (cm)	Order Code with Guide Wire	Guide Wire Length (cm)	Guide Wire Diameter	Guide Wire Tip	Guide Wire Stiffness
07	12	26	TU-201226-D	125	0.035"	straight	standard
07	12	28	TU-201228-D	125	0.035"	straight	standard
07	12	30	TU-201230-D	125	0.035"	straight	standard

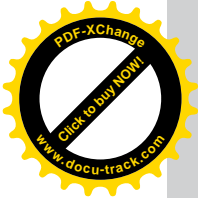

URETERAL STENT-SET for supravvesical urine drainage

SET COMPONENTS



Ureteral Stent

- Polyurethane (PUR)
- 8 drainage openings at distal end for constant urine flow
- Open tip
- Markings for exact positioning
- Radiopaque
- Application e.g. after neobladder reconstruction surgery



Adapter

- Cannula adapter for drainage or attachment of the urinary drainage bag (Luer Connector)

Open tip					
Size (FR/CH)	Length (cm)	Material-Typ		Order Code without Guide Wire	Max. Guide Wire Diameter
4.7	70	White-Star	○	WEB-120570-NEO	0.028"
06	70	Green-Star	●	WEB-120670-NEO	0.035"
08	70	Green-Star	●	WEB-120870-NEO	0.035"